

# Attachments and Pre-fabricated Castable Components

CATALOG/TECHNICAL MANUAL


for Dentists and Dental Technicians



13<sup>th</sup>Edition

World Leader in Dental Attachments

### **DENTIST COURSES AND UNIVERSITY PROGRAMS**

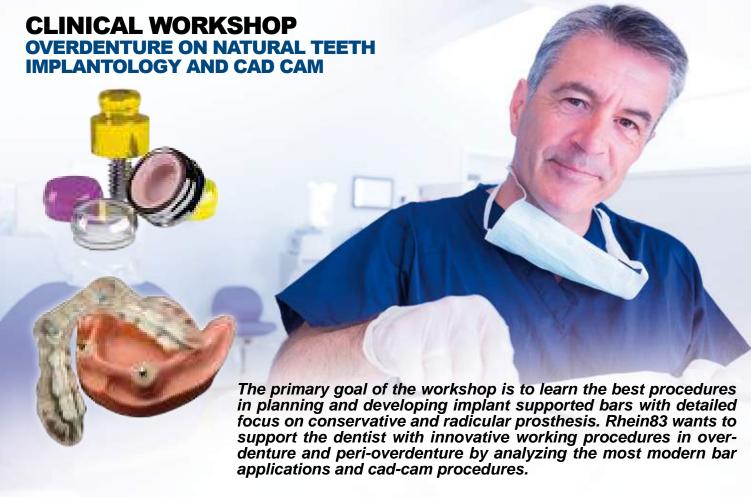

### **UNIVERSITY PROGRAMS**

# REMOVABLE PROSTHESIS MASTER COURSES, TRADITION AND INNOVATION OF THE RETENTIVE SYSTEMS









Courses dedicated to universities presenting innovative solutions and procedures in planning the prosthetic projects. Functional, aesthetical phonetical evaluation of the patient by considering the social conditions and background. Real clinical cases presentation and analysis supported by live working procedures on models with students from universities worldwide. Cultural interchange programs with international universities, post graduate degree programs, international contests and much more!

### **DENTAL TECHNICIAN COURSES**

### **BASIC LEVEL**

Introduction to the Rhein83 techniques in intra-coronal and extra-coronal prosthesis. Innovative procedures allowing to reduce working times and costs by using pre fabricated castable components. Direct overdenture concepts in implantology on all implant brands and platforms.







### **MASTER LEVEL**

Deeper insight into the themes presented during the basic course with special focus on implant prosthesis and new digital cad cam working procedures. Simple and useful solutions in complex implantology clinical cases.









### **RHEIN83 BIRTH, GROWTH AND EVOLUTION**


Metallic spherical attachments exist since many years. But these attachments were not widely accepted, by the dental professionists. Then came the idea to render these mechanisms elastic! A smoothed head and the elastic cap are the result of these innovative changes; today this technique is amongst the most widely used. Rhein83 has been in business since 1983 and today these products have been copied throughout the entire world, copies that in many cases reflect the forms of the objects but not the materials they are made from, and therefore it significantly changes the functional result. Research is not only oriented towards the study of new products, but also continually trying to perfect those that have been used for many years. Dental attachments are small mechanisms subjected to continuous movement, stresses and oral changing, requiring periodic maintenance and revisions. Some products in this have been made for maintaining and restoring the functionality, to all the prostheses, directly while they are in the mouth of the patients. The commitment of Rhein83 with its knowledge and skills continually being enriched by the contributions of dentists and laboratory technicians, is to be able to improve the actual standards and develop new products by means of original projects.

Ezio Nardi



#### **RESEARCH AND INNOVATION TODAY**

By over 34 years Rhein83 is continuously innovating the dental attachments world with materials and designs allowing to satisfy the technical requests of the dental specialists.

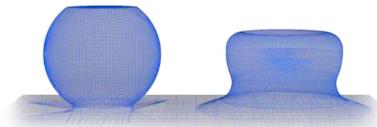

#### **NEW OT EQUATOR PROFILE**

Evolution from the sphere to the semi-sphere, reduced dimensions allowing the same stability and functionality!

#### TECHNICAL INNOVATIONS AVAILABLE TO ALL!








# **GENERAL INDEX**

RHEIN83 BIRTH, GROWTH AND EVOLUTION	2
GENERAL INDEX	3
FRICTIONS AND RETENTIONS CONCEPT	4
FEMALE CAPS ASSORTMENTS	5
OT EQUATOR CASTABLE	6-7
OT EQUATOR FOR IMPLANTS AND SMARTBOX	8-9
OT EQUATOR ELASTIC SEEGER	10-11
OT CAP SINGLE THREADED SPHERES	12-13
OT CAP & OT CAP TECNO - COMBINED PROSTHESES	14-15
OT BOX MONO	16
OT STRATEGY - COMBINED PROSTHESES	18-19
OT STRATEGY/STEADY	20
OT STRATEGY & OT CAP PROSTHETIC PROJECT	21
SINGLE SPHERES - OT CAP CASTABLE - OT CAP TITANIUM + TIN DIRECT SYSTEM OVERDENTURES	22-23
S.P.L. TITANIUM POSTS FLEX - BLOCK DIRECT SYSTEM OVERDENTURES	24-25
OT BOX, CLASSIC - SPECIAL - CAST REINFORCEMENTS WITHOUT MODEL DUPLICATION	26-27
OT REVERSE 3 DIRECT SYSTEM OVERDENTURES	28-29
RECONSTRUCTIVE SPHERES: CONCAVE SPHERE - OT EQUATOR	30
RECONSTRUCTIVE SPHERES: SOLID SPHERE	
OT BAR MULTIUSE	32-33
OT VERTICAL	34-35
OT UNILATERAL	36-37
OT LOCK LOCKING PIN	38-39
IMPLANT OVERDENTURE ATTACHMENTS: SPHERO FLEX - BLOCK, DIRECTIONAL RINGS	40-41
IMPLANT OVERDENTURE ATTACHMENTS: UNIVERSAL "ANTI-UNSCREWING" SYSTEMS	42
MINI PARALLELOMETER DEVICE WITH MODEL HOLDER BASE AND CUFF HEIGHT MEASURER	43
IMPLANTOLOGY: BROKEN SCREW EXTRACTOR FOR IMPLANTS FOR REMOVAL OF BROKEN IMPLANT SCREWS	44-45
POLISH BURS KIT BY CARLO BORROMEO	46
INSTRUCTION AND TECHNICAL ADVICE	47
ACRYLIC DEMONSTRATION MODELS	48
PRODUCT RANGE - SIZES AND DIMENSIONS	49-50
KITS AND CODES	4-55-56
RHEIN83 WORLD WIDE	57



### **COMPARISON OF RIGID CAPS vs. ELASTIC CAPS**

#### Characteristics and retentive functionality

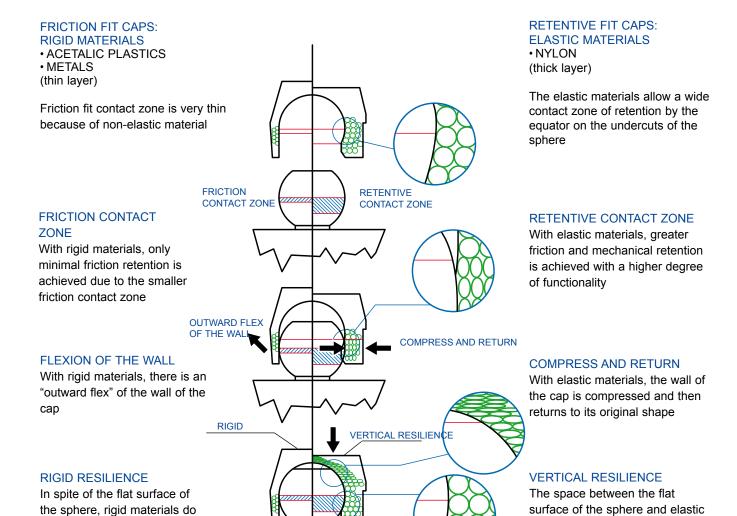

#### **RHEIN83 - DESIGN AND FUNCTION**

not allow vertical resiliency

Rhein83 continues to manufacture female caps with elastic retention with the intention of eliminating as much vertical stress and trauma to the restoration as possible. For Rhein83 the important thing is to make a system of components available to the dental technician and dentist that will allow for the fabrication of a rigid, shock absorbing or resilient prosthesis. With the use of elastic retention, the function of Rhein83 attachments are extended.

With overdenture prosthetic devices or cases involving edentulous saddles, resiliency can be controlled with a wide range of retentive caps that have various levels of elasticity and retention.



cap allows for vertical resiliency

and reduces stress

## **CLASSIC CAPS** SIZES AVAILABLE: NORMAL AND MICRO Retentive cap colors and retention

#### **CLEAR CAPS** STANDARD RETENTION





#### Slightly Elastic

Maximum suggested time of duration in mouth: 12 months Retention in grams: Normal 1300g / Micro 1100g

#### PINK CAPS SOFT RETENTION



Maximum suggested time of duration in mouth: 12 months Retention in grams: Normal from 900g / Micro 800g

#### YELLOW CAPS

**EXTRA SOFT RETENTION** 



#### Verv elastic

Maximum suggested time of duration in mouth: 12 months Retention in grams: Normal 500g / Micro 450g

#### **GREEN CAPS ELASTIC AND GUMMY**



**Characteristics** Maximum suggested time of duration in mouth: 12 months Extremely elastic retention, "GUMMY" type. Minimally hydroscopic, with a good adhesion on the sphere.

Retention in grams: Normal 350g / Micro 200g

#### **FXTRA RESILIENT** GOLD CAPS SLIGHTLY ELASTIC



Characteristics Maximum suggested time of duration in mouth: 12 months To be used in overdenture prostheses, where resilience and vertical movements are necessary.

Retention in grams: Normal 500g / Micro 450g

**EXTRA RESILIENT** SILVER CAPS **ELASTIC AND GUMMY** 



Characteristics Maximum suggested time of duration in mouth: 12 months To be used in overdenture prostheses, where a vertical movement is necessary and a light initial retention is requested. Retention in grams: Normal 350g / Micro 200g

#### PROCESSING CAPS





#### Characteristics

Caps to be used only for laboratory processing.

#### **TITAN CAPS** NYLON CAPS WITH INTERNAL TITANIUM RING





Characteristics Maximum suggested time of duration in mouth: 12 months Extremely durable. To be used especially in combination with pre-fabricated spheres such as titanium spheres, concave spheres, etc. Retention in grams: Normal 1500g / Micro 1300g

#### **UNDERSIZED INTERNAL** DIAMETER CAPS STANDARD RETENTION





Characteristics Maximum suggested time of duration in mouth: 12 months Internal diameter reduced (Normal 2.2mm | Micro 1.6mm), for 2.25mm - 1.6 spheres

Retention in grams: Normal 1300g / Micro 1100g

#### UNDERSIZED INTERNAL DIAMETER CAPS



Characteristics Maximum suggested time of duration in mouth: 12 months Internal diameter reduced (Normal 2.2mm), for 2.25mm spheres


Retention in grams: Normal 900g

### SOFT RETENTION **UNDERSIZED INTERNAL**



**Characteristics** Maximum suggested time of duration in mouth: 12 months

**DIAMETER CAPS EXTRA SOFT RETENTION** 

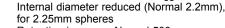

**UNDERSIZED INTERNAL** DIAMETER CAPS

Retention in grams: Normal 500g

**ELASTIC AND GUMMY** 

Characteristics Maximum suggested time of duration in mouth: 12 months Internal diameter reduced (Normal 2.2mm | Micro 1.6mm), for 2.25mm - 1.6 spheres

Retention in grams: Normal 350g / Micro 200g

#### STAINLESS STEEL AND TITANIUM HOUSING FOR CAPS, PRE-FABRICATED, NORMAL AND MICRO SIZES



The new stainless steel housing design offer reduced size and additional stability, it can be enbodied directly in the resin, welded or bonded to the frame. The new design is also available in titanium.

# **EXTRA RESILIENCY FUNCTIONALITY** Extra resilient caps, normo and micro size, will allow to absorb elevate masticatory forces without creating any damage to the implant or root.

### **OT EQUATOR CASTABLE**

# Single Attachment for Overdentures





RETENTIVE CAPS

OT EQUATOR

OT EQUATOR CASTABLE MALE ATTACHMENT



**IMPRESSION** TRANSFER pick-up impression



**IMPRESSION** TRANSFER individual tray



TITANIUM HOUSING



VIOLET CAP RIGID RETENTION (2.7Kg)

STAINLESS STEEL HOUSING



CLEAR CAP STANDARD RETENTION (1.8Kg)



PINK CAP SOFT RETENTION (1.2Kg)



YELLOW CAP EXTRA-SOFT RETENTION (0.6Kg)



BLACK CAP **PROCESSING** 



STAINLESS STEEL ANALOG FOR PLASTER MODEL






PARALLEI OMETER MANDREL


OT EQUATOR CAPS







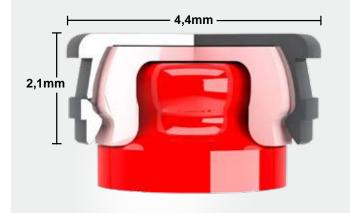

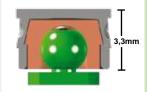



If additional retention is needed to secure the prosthesis, OT Cap Normal retentive caps and metal housings can be placed over any OT Equator Profile spheres. The prosthesis will be retained in the same way and the connection will be more rigid. Only the dimension of the attachment will be changed.



SEVERE DIVERGENCY MAY REQUI-RE THE OT EQUATOR IN COMBINA-TION WITH A CASTABLE UCLA

# LABORATORY

#### OT EQUATOR CASTABLE = INDIRECT TECHNIQUE



Use separating material on the stone model in the prepared areas to receive the castable posts.



Use longer castable posts in the root channels for easy removal. Reline with castable resin, for higher accuracy.



Place posts and finish margins with composite material. Once resin is cured, cut posts to the required length at root level.



Position OT Equator on the occlusal surface with the paralleling key and continue waxing technique.



OT Equator in the final position. The waxup has been completed.



For the best results, create the casting with an alloy that has a vickers hardness of 220 or greater.

#### BUILD UP THE FRAME DIRECTLY ON MASTER MODEL



The plaster model with the OT Equator analog in position. The stainless steel housing and black processing cap are also visible.



Apply a thin layer (.5mm) of wax to the model. Fill the undercuts on the stainless steel housing and attach the connectors.



Connect the parts using a castable resin. Be sure to cover the stainless steel housing.



Add sprues to the framework and remove it from the model. Be sure that the stainless steel housing does not remain inside. The framework is now ready to be invested.



Cast the metal frame and verify the position on the model.



Use composite to bond the stainless steel housing to the frame.




The metal frame with the stainless steel housing in place.



The finished prosthesis on metal frame. After processing, the black caps are replaced with pink caps.

### **OT EQUATOR FOR IMPLANTS**


### Low Profile Titanium Abutment

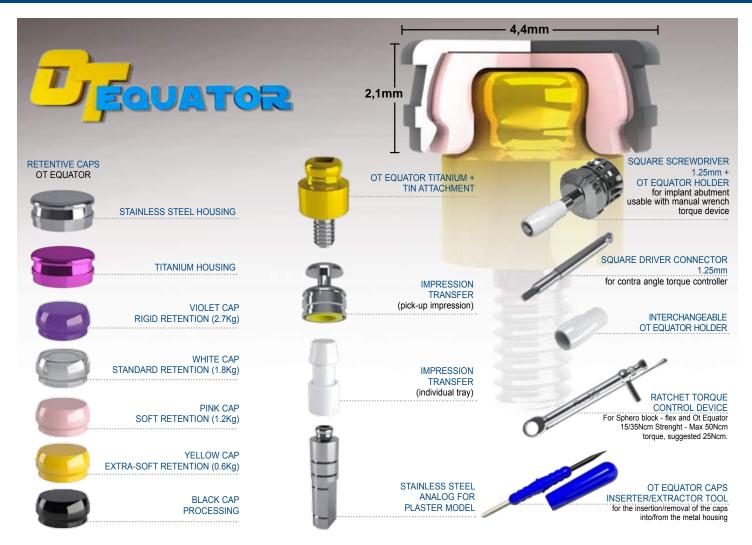

The unique design and exceptionally low 2.1mm profile of the OT Equator 4 in 1 System provides exceptional stability and superior retention when compared with other attachment systems. Due to its lower radius, OT Equator is indicated to correct divergence up to 25 degrees between implants without affecting the functionally of the elastic nylon cap. Caps are available in a wide variety of retention levels. ATTENTION; Where implant divergence exceed the maximum 25 degrees, Sphero Block and Sphero Flex are recommended case plan options. See Sphero Block and Sphero Flex page 40-41



# the **self-aligning** Ot Equator Housing

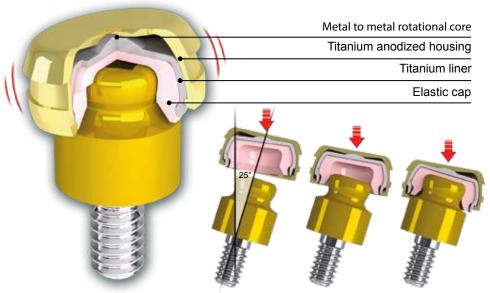



TITANIUM HOUSING WITH BLACK CAP



SMARTBOX BLACK POSITIONING CAP

Passive insertion reduces trauma

Correct divergency up to 50°



#### ATTACHING THE CAPS IN CLINIC



Select the OT Equator with the appropriate cuff height. Screw the OT Equator into the implant.



Place the protective disk over the OT Equator. Then, place the stainless steel housing with cap on the attachment.



Verify the positioning of the prosthesis before bonding the stainless steel housing.



On the prosthesis, fill the implant sites with a self curing resin and insert into the patient's mouth.



Remove the prosthesis and verify that the positions of the attachments are correct.



Remove the protective disks.



Carefully trim away the excess The completed prosthesis. resin.



#### **IMPRESSION TRANSFER**



Place the impression coping on the OT Equator.



analog Insert the into impression coping and pour the master model.



Add sprues to the framework and remove it from the model. Be sure that the stainless steel housing does not remain inside.

**BUILD UP THE** FRAME DIRECTLY (for the full technique go to page 7)



The metal frame with stainless steel housings bonded in place.

#### CHAIRSIDE PROCEDURE FOR SMARTBOX POSITIONING



appropriate cuff height. Screw the OT Equator. OT Equator into the implant.



Select the OT Equator with the Position the protective disk over the



Fully engage SMARTBOX with Black Fill the space corresponding to cap securely onto OT Equator.



the housings with self curing resin. Insert the prosthesis into the final position.



Once the resin has cured, remove the protective disk.



Remove excess resin with bur and polish for passive connection.



Remove SMARTBOX black cap with cap extractor tool.



Using the cap insertion tool, select 1 of 4 Ot Equator femal caps for desired retention.



### **ELASTIC SEEGER**

### Passive bar connection













The purpose of the OT Equator "seeger" system is to create a passive connection for implant supported bars. The elastic seeger will correct small imperfections created by the chairside impression technique or laboratory casting process. This reduces the risk of the implant bar to not seat passively.

#### POSITIONING SYSTEM WITH BAR "ELASTIC SEEGER"



OT Equator titanium attachments screwed into the implants. The elastic seeger system will be used to position the bar.



The cast bar in position. Insert the PEEK elastic seeger ring into the cylindrical space.



Using the insertion tool, push down the PEEK elastic seeger ring until it is fully seated.



PEEK seeger ring in position, ready to screw the titanium locking screw.



After the elastic seeger ring has been inserted, lock the bar into place using the titanium locking screw, (Torque suggested 15 Ncm)



The finished bar secured in the mouth. passive connection has achieved due to the elastic PEEK seeger rings.



The completed prosthesis. For best results a reinforced superstructure is always recommended.



In case of a future check, the special internal design of the PEEK seeger ring allow the self extraction together with the titanium locking screw

### WAX-UP OF THE BAR DIRECTLY ON MODEL MASTER



Screw the OT Equator attachments into the implant analogs.



Position the seeger castable cylinders, followed by the red plastic seeger for laboratory use on the attachments (Thinner part lower). Screw the titanium sealing lid into position. Do not overtichten not overtighten.



OT EQUATOR castable attachments or equation attachments are placed on the connecting bar creating a "balance" with the removable prosthesis. Alloys with a Vickers Hardness of 240 or greater are recommended for carting recommended for casting.



Connect the castable abutments with



The cast bar in position on the model.



The cast framework in position. Undercuts on the stainless steel housing can be blocked out using composite material to maintain a passive connection.



Fit and stability of the prosthesis can be regulated using nylon caps. Various levels of retention are available.



The final prosthesis.

#### OT EQUATOR SEEGER



Insertion of the Ot Equator seeger inside the self-extracting peek ring already in position.



Like the titanium seeger screw, the Ot Equator seeger also removes self-extracting peek ring during unscrewing.



2 threaded Ot Equator for cam and 2 Ot Equator seeger in position, in case of divergence (the Ot Equator seegers follow the direction of the implants) the use of Smart box housings is recommended.



## **INTERCHANGEABLE THREADED ATTACHMENTS**

# with threaded sleeve system





#### THREADED SLEEVE FOR BONDING



OT CAP SLEEVE **SPACERS** Normal/Micro

**SCREWDRIVER** HFX 0.9 mm



**SCREWDRIVER** HFX 13 mm

#### RETENTIVE CAPS OT CAP



Standard













Black Processing



OT EQUATOR SQUARE HEAD



thread 1.6 mm



#### THREADED SLEEVE FOR BONDING



SCREWDRIVER OT EQUATOR SQUARE HOLDER



S.STEEL OR **TITANIUM HOUSINGS** OT EQUATOR



OT EQUATOR SLEEVE SPACER

#### RETENTIVE CAPS OT EQUATOR











Black





















Processing









#### STEP BY STEP THREADED SLEEVE BONDING PROCEDURE



Once the bar has been connected with wax, create an area where the attachment spacer will be placed.



Apply separator to the base of the attachment spacer and postion using the parallelometer key.



With the attachment spacer in position, complete the wax-up design.



Carefully remove the attachment spacers and proceed with the normal casting procedure.



Screw the threaded attachment of choice (Micro Ball shown) into the threaded (Micro sleeve.



Place the assembled attachment into the parallelometer key. Use a self curing metal to metal bonding composite on the sleeve and in the cylinder.



After the composite is cured, remove any excess material.




Unscrew the attachment to verify if the threaded sleeve is securely bonded in place.



The finished bar complete with attachments.

### **3 ATTACHMENT OPTIONS**

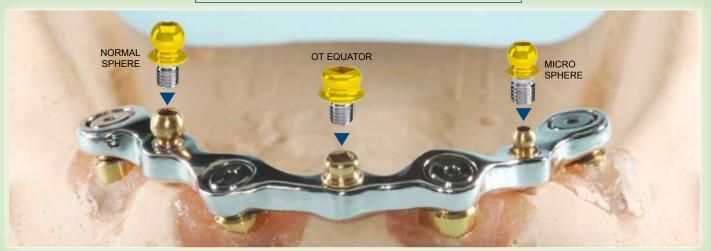

THE TECHNIQUE IS THE SAME FOR ALL THREE OPTIONS

## **EXTRACORONAL CASTABLE ATTACHMENTS** OT CAP - OT CAP TECNO



OT Cap is a resilient distal extension attachment. It is indicated to be used with combined prostheses and removable partial dentures.

For treatment plans that require a rigid substructure with milling and adequate counter attachments, OT Cap functions as a stabilizing retentive connector. In addition, for treatment plans which require resiliency, OT Cap provides a "Cushion Effect" similar to a shock absorber. This is achieved by the design of the sphere in conjunction with the elastic retentive caps.

The OT Cap Tecno consists of a titanium sphere and ring that is incorporated into the nylon cap which has been machined with a tolerance that assures high precision. While fabricating the prosthesis, the Tecno titanium sphere is not exposed to any of the risks associated with the laboratory fabrication procedures and ceramic firing cycles.





OT EQUATOR CAPS

NSERTER/EXTRACTOR TOOL e insertion/removal of the caps into/from the metal housing

> View of the Ot Techno system, Normo or Micro sphere can be used with the same threaded sleeve.

#### OT MONO BOX



OT BOX MONO: The positioning ring to be inserted on the sphere before model duplication.

REGULATING TOOL FOR

RETENTION

Normal / Micro

# **COMBINED PROSTHESES**

### with extracoronal castable attachments





















#### OT CAP CASTABLE



Cut the plastic bar and use only the section that you need.



Using the mandrel, position the spheres in parallel. Complete the wax-up with a "ledge" along the crown. The "ledge" must not be lower than the sphere.



cast crowns. It suggested to use a retentive cap to protect the sphere from any damage.



The cast attachment. The "ledge" along the crown helps select and redirect the vertical loads.



Using the mandrel, position the Ot Tecno castable extension in parallel, Complete the wax-up with a "ledge" along the crown and cast.



Place the assembled attachment into the parallelometer key. Use a self curing metal to metal bonding composite on the sleeve and in the cylinder.



After the composite is cured, remove any excess material.



Unscrew the attachment to verify the threaded sleeve is securely bonded in place.

#### CAST HOUSING WITH DUPLICATED MODELS



The OT Cap positioning ring on the sphere.



The duplicated model in investment.



The OT Mono Box castable housing in position and housing in position and incorporated into the final wax design.




The final OT Mono Box casting with retentive caps inserted into the housing.

The castable OT MONO BOX reproduces the shape of the housing which incorporates the retentive cap into the framework. Use the OT CAP insertion tool to place the retentive cap into the housing.



# **CASTABLE** Customized solution for frames with single castable sphere housing for caps

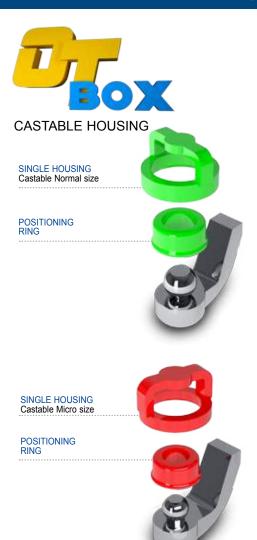

#### **HOUSINGS:**

#### STAINLESS STEEL - TITANIUM

The new stainless steel housing design offer reduced size and additional stability, it can be embodied directly in the resin, welded or bonded to the frame.

The new design is also available in titanium.



When vertical space is limited, use reinforced pins to reduce the risk of breakage of the denture teeth.

#### **SOLUTION A**



Place a piece of .5mm calibrated wax over the wax-up design for additional protection.



The finished casting with retentive cap

#### SOLUTION B



pins are added reinforcement of the denture acrylic as well as additional retention for the denture teeth.



The final cast housing with reinforced metal pins.

#### STAINLESS STEEL PRE-FABRICATED HOUSINGS

For bonding or soldering to the frame

To obtain the right position use the POSITIONING RINGS. NORMAL and MICRO sizes are available.











### **CERTIFICATIONS**

Rhein83 continues to be the world leader in spherical attachments and implant components. Largely due to continuous research and development, active participation in exhibitions as well as providing practical hands-on technical training for dentists and dental laboratory technicians. In addition, the company utilizes state of the art technology to constantly develop new products and improve existing product désign as well as promote product awareness.

Rhein83 attachment systems are technically supported in over 75 countries worldwide.













#### **CERTIFICATIONS**

Rhein83 started a "Quality System" path in 1994, obtaining ISO 9001 certification since 1996. In the same year, the company obtained the first product accreditations from the FDA/USA and continued the path with the certifications: 1998-93 / 42 EEC - ISO 9001: 2015 - ISO 13485: 2016. Other conformities and registrations have been obtained in many foreign countries, including: Russia, Canada, Brazil, Korea, etc. The "regulatory" activity is always monitored and constantly updated, for the maintenance of current certifications, their renewal and adaptation to the new rules for medical devices.



Claudia Nardi Administrator



Gianni Storni Administrator and R&D director



### **CASTABLE VERTICAL ATTACHMENT MICRO**











#### CAPS



Clear • Standard



Pink • Soft



Yellow • Extra Soft



INSERTER/ **EXTRACTOR TOOL** 

MANDREL

OT STRATEGY CAPS

**PARALLELOMETER** 

Black • Processing











CAPS

















cap housing on the frame



the OT Insert Strategy male into the mandrel and place position with base of attachment in contact with the stone.



The entire cap must be covered with a thin layer of wax during the frame wax-up procedure.

#### PARALLELOMETER KEY PROFILE





SIDE A: For SPHERE positioning SIDE B: For STEADY positioning

#### REINFORCEMENT FOR THE SPHERE







- Increased shear force strength
- Prevents rotation of female cap
- Increased lateral stability

OT Strategy from Rhein83 is a vertical micro-sized 1.8 mm castable sphere that is placed distally on abutments for removable partials or utilized in implant bar combination case design. The male component is designed with an additional support strut located under the sphere, increasing strength and preventing rotation of the female cap during paralleling. The optional Steady, when connected to OT Strategy, provides lateral stability without any additional milling

OT Strategy caps are available for both duplication and fabrication using a stainless-steel housing technique. Rhein83 caps are manufactured from an elastic material that increases the contact zone with the sphere, giving mechanical and friction retention. Caps are color-coded indicating five levels of retention. Tools for paralleling, inserting, and removing caps are available.



Once the casting complete, proceed to use the cap and the prefabricated **STAINLESS** STEEL HOUSING. The housing can bonded or laser welded to the frame. In addition, it can also be used for direct chairside procedures.



For best results during **DUPLICATION** the TECHNIQUE, it is suggested to use the YELLOW retentive cap.



# **MIXED PROSTHESES**















#### **DUPLICATION TECHNIQUE: USING CASTABLE HOUSING**



OT Strategy casting is complete with mandatory lingual milling to accept partial bracing arm. accept



Yellow retentive cap is placed on the sphere and the model is ready for duplication. Use wax to remove any undercuts.



Model is duplicated and the shape of the cap is reproduced.





Insert the black cap into the skeletal cast frame cast partial with the OT Strategy Insertion Tool.



Frame is complete placed on the model.



Using the insertion tool, insert the cap.



The finished prosthesis.

#### WELDING TECHNIQUE: USING PRE-FABRICATED STAINLESS STEEL HOUSING



Crown and OT Strategy attachment cast. Positioning ring and housing.



Wax-up on the duplicated model.



Positioning ring on the sphere.



First Option: Stainless Steel Housing welded to frame.



Stainless Steel Housing in position on the attachment.



Second Option: Stainless Steel Housing bonded to frame with anaerobic selfcuring resin.



ATTENTION: Insertion of the cap from the mesial.

# **CASTABLE VERTICAL** MICRO ATTACHMENT STRATEGY + OPTIONAL STEADY












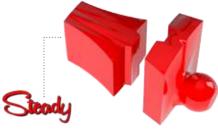

#### TECHNIQUE WITH STANDARD BASE









Lute the two parts together using an adhesive and insert the sphere into the mandrel of the parallelometer.

The Steady can be used with it's original height or it can be shortened and modified to accommodate the adjacent tooth and ridge.

Finish the wax-up give the Steady necessary shape for duplication in the sphere.





The frame wax-up.



Steady + standard base

**Optional = STEADY** 

TECHNIQUE WITH LONG BASE



The finished casting.

Steady + long base

Lute the Steady to the Long base. Be sure to position the two parts according to the resorption of the ridge. Position the attachment as close to the ridge as possible. Fill the space between the Steady and the ridge with wax

The finished attachment design. The Steady has been adapted to the contour of the ridge.

The castable Steady is an optional conical shaped support intended for use in cases where milling is not performed. Steady can be used with the OT Strategy Standard or Long base.







It is an object in line with the philosophy of the personalization of each single prosthesis and

Crown and Steady for duplication and retentive cap on the sphere. Cast framework seated on the model.


Finished prosthesis.

with both the OT Strategy bases; Standard or Long and offer various technical solutions.

When the STEADY base is utilized it provides superior lateral support when milling is not indicated. In the case of free saddles, the STEADY base avoids movement in all directions during mastication.

# **OT STRATEGY & OT CAP**

# Case design options

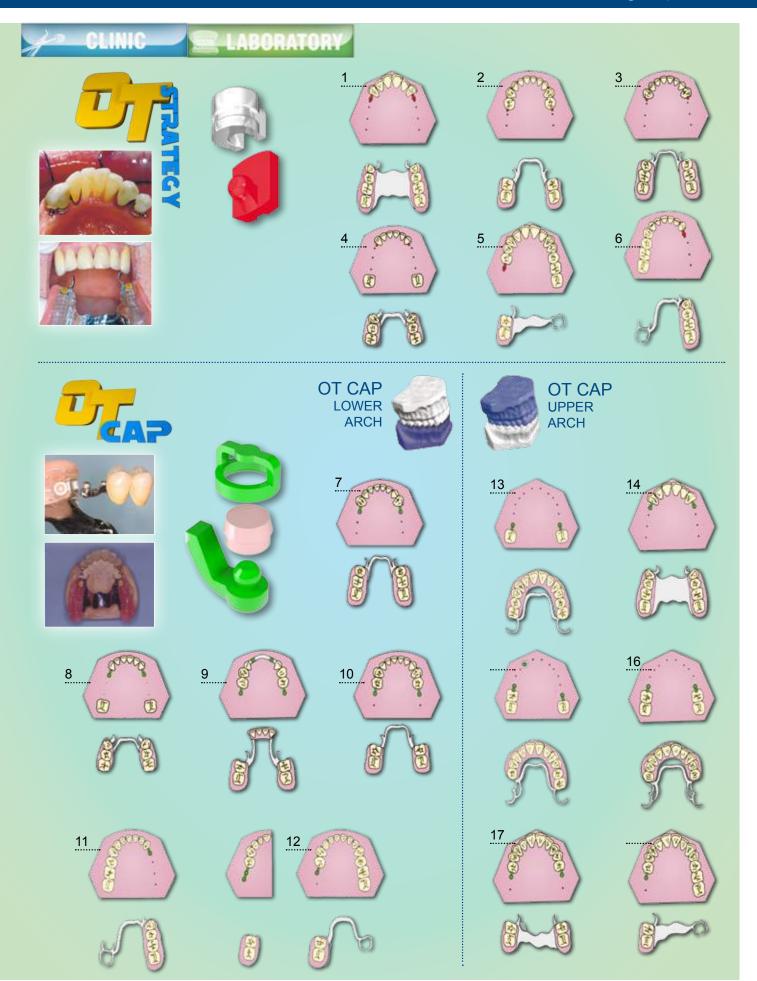

### **SINGLE SPHERES OT CAP**



#### STAINLESS STEEL AND TITANIUM HOUSINGS

Normal / Micro size for curing welding or bonding





**ELASTIC RETENTIVE CAPS** 

Normal / Micro

Clear • Standard

Yellow • Extra Soft


Green • Elastic

Pink • Soft



SINGLE SPHERES TITANIUM + TIN

1600 Vickers Hard FOR WELDING OR BONDING

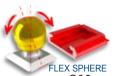





impression copings and pour FIXED SPHERE the model. NORMAL SIZE



TRANSFER IMPRESSION TECHNIQUE

Different

Stone model with analogs in

#### **CASTABLE SINGLE SPHERES**



NORMAI Green Ø 2.5 mm



Ø 1.8 mm

Black • Processing

Undersized caps for worn or damaged spheres are also available. They are also compatible with 1.7mm and 2.2mm spheres

See parts list for item codes and descriptions.

for impression of the root canals



FIXED SPHERE

Ø 2.5 mm





#### ATTENTION:

These attachments can be cast with all types of alloys, but it is important to use a metal with a high Vikers hardness in order to avoid the risk of wearing the spheres.



**IMPRESSION** COPING Normal / Micro

MOOSER BURS



Normal











OT EQUATOR CAPS INSERTER/EXTRACTOR TOOL for the insertion/removal of the caps into/from the metal housing

The design of the sphere with a FLAT head in addition to the spherical inner surface of the elastic cap, permits vertical movement during mastication. Rhein83 female caps are manufactured out of a special nylon material that remains stable and continues to function in the oral cavity over long periods

Clinical data is available showing that stability is obtained with a minimal amount of trauma.





## **OVERDENTURE PROSTHESIS**

## **Indirect System**


#### IMPRESSION OF ROOT CANALS



Prepare the roots.




Apply adhesive to the post.

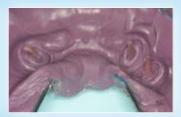

Impression with elastomer.

#### ATTENTION:

To obtain proper function, it is important to mill the resin with a bur to create a space (highlighted in blue) between the root and the prosthesis.

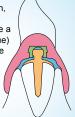

#### OT CAP - EMBODING STAINLESS STEEL HOUSING TO DENTURE



Protective discs on the cast metal spheres.



Fill the space corresponding to the housings with self curing resin. Insert the prosthesis into the final position.



Once the resin has cured, remove the disc and trim the excess material around the housing.



Finished prosthesis.

#### OT CAP - CASTABLE SINGLE SPHERE TECHNIQUE



Insert the castable plastic post into the prepared root cavity.



Cut the post to the level of the root and remove the sphere.



Position the single spheres in parallel with each other.



Cast post and sphere. It is also possible to place the sphere off center in respect to the long axis of the post.

#### OT CAP - TITANIUM SINGLE SPHERES + TIN FOR CURING WELDING OR BONDING



Wax-up the root cap. Insert the titanium sphere into sliding base and position it on the root cap.



Wax-up with titanium sphere in position. Do not cover the "open" side of the base with wax.



Remove the titanium sphere from the base before attaching sprue.



The finished wax-up with sprue. The root cap and post is ready to be invested



Using the tool, check the fit of the cast cap by inserting the sphere into the base.



Titanium sphere inserted in the cast root cap base.



Bond the titanium sphere to the base using anaerobic or self curing composite material.



Finished root cap. The sphere is bonded and locked in position.

### **PIVOTS FOR DIRECT OVERDENTURE**



**PIVOT FLEX TITANIUM +TIN** 1600 Vickers Hard "self-paralleling" sphere



STAINLESS STEEL AND TITANIUM HOUSINGS Normal / Micro size curing welding or bonding











#### SUPER-RESILIENT CAPS





**TITANIUM** PIVOT BLOCK

(Ø 2,5 mm/Ø 1,8 mm)









LOW PROFILE **TITANIUM + TIN PIVOT** 

(over 1600 Vickers)







TITANIUM HOUSING



SMARTBOX HOUSING



VIOLET CAP RIGID (2.5Kg)



TRANSPARENT CAP STANDARD (1.8Kg)















PIVOT BLOCK TITANIUM WITH STATIONARY SPHERE NORMAL Sphere







OT EQUATOR PIVOT LOW PROFILE TITANIUM PIVOT 3 lengths: 10,9,7mm







The Pivot Flex line of titanium posts was developed as an economical solution for direct "in root" supported overdentures. The self-aligning Pivot Flex post features a rotating ball with a 2.5 mm diameter and is indicated for divergent roots. When the posts are used with directional rings to align retentive caps before the resin curing stage, the insertion of the denture is easy and trauma-free. The Pivot Block line of milled titanium posts has a stationary ball and can be used for a temporary or as a permanent solution. The Pivot Block titanium posts are available in 2.5 mm and 1.8 mm sphere diameters. The Rhein83 elastic caps ensure optimal retention and function while minimizing wear. There are five levels of retentive caps, including extra resilient caps for precarious root situations. The levels of retention are identified by different colored caps.

#### **DIRECTIONAL RINGS - FOR FIXED AND ROTATING SPHERES**



Pivot Flex posts in divergent roots



Nylon caps without directional rings. Caps are not supported in the same horizontal plane. without



Nylon caps with directional rings. Caps are now supported in horizontal plane. in the same

### DIRECTIONAL RINGS



14°

#### PIVOT BLOCK - FOR TEMPORARY OR PERMANENT ECONOMICAL SOLUTIONS



Pivot Block cemented with oxyphosphate cement for a temporary solution.



To remove the post from the root, grasp the sphere with the pliers and rotate carefully in both directions.



Due to the conical shape and smooth surface, the post is removed easily.



For permanent solutions. create notches in the post and roughen the surface before cementation.







#### TITANIUM PIVOT BLOCK AND OT EQUATOR: PERMANENT FIXATION IN THE PATIENT'S MOUTH



Prepare the root by the mucosal level and adjust the radicular cavity by using a Mooser Bur with the proper dimensions



Alginate impression: attachment placements in evidence.



Fill-up the radicular cavities with proper composite cements, insert than the spherical titanium pivots.



Place the protective disks between the directional rings and the retentive caps. Feel with self curing resin and than place the prosthesis in the patient's mouth.



Cemented micro block pivot in position, retentive notches were applied to support the permanents fixation.



When the resin will be hard enough remove the protective disk and clean up any excess of resin.



Place the directional rings in position between the roots and retentive caps. Proceed by taking the imprint.



Completed prosthesis.

### **CASTABLE BAR CAP HOUSINGS**

Clear • Standard Pink • Soft

Yellow • Extra Soft

Green • Elastic Black • Processing

Gold • Slightly Elastic

**IMPRESSION** 

**COPING** Normal / Micro

Silver • Elastic and Gummy



Universal castable bar for joining the OT BOX housings







SPECIAL BARS NORMAL - Green MICRO - Red







**LARGE O** 2 **NORMAL + CONNECTOR** 




MICRO + CONNECTOR



Cast reinforcement wax-up on the master model without duplication.

The OT Box Large casting compensates for the distance between the cap and the housing. It is manufactured to reposition the cap charside into the



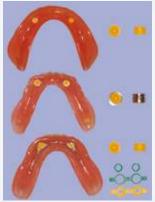

A fracture is more likely to occur where the overdenture attachments are inserted in a prosthesis fabricated entirely of resin. With a cast superstructure reinforcement, the denture will be less likely to fracture. Fast and simple, the OT Box bar components are used to fabricate the superstructure directly on the master model, eliminating duplication and saving time. A nonprecious or chrome cobalt alloy is recommended for . best results.

It is recommended that all nylon caps are inserted into a stainless steel housing or cast reinforced frame. The stainless steel housing offers a considerable advantage when the cap has to be removed and replaced for routine maintenance or repositioned. Adjustments or repairs can be performed chairside quickly and easily.

Option 1: OT CAP OT Cap cured directly into the prosthesis.

Option 2: OT Cap + Stainless Steel Housing OT Cap with housing cured directly into the prosthesis or bonded into frame.

Option 3: OT Cap + OT Box OT Cap inserted into OT Box cast reinforced frame.

#### LABORATORY REQUIREMENTS FOR THE MASTER MODEL

When a new denture is being fabricated utilizing existing spheres, the dentist must provide the laboratory with an impression using the YELLOW CAP. The laboratory will place the analog into the cap and pour the stone model.





## **CAST REINFORCEMENT IN ACRYLIC DENTURES**

# without duplication of the model

#### IMPRESSION WITH POSTS FIXED IN THE MOUTH



Titanium posts cemented into the root.



Before taking the imprint place the transfert over the spheres supported by the proper directional ring.



Insert analogs into the impression copings and pour the model.



Stone model with analogs in place.



Plaster model with metal-fused components.





#### DIRECT WAX-UP ON THE MASTER MODEL



OT Box Classic. Glue the two OT Box bars together.



Separate the housing from the OT



"ONE-PIECE" MONO **BAR** OT BOX SPECIAL is a "one-piece" mono bar. Separate the bar and use only the section needed.






Apply a layer of wax on the ridge. Create three holes in contact with the stone model. Place the positioning rings over the spheres. Be sure to place the ring with the "flared" end towards the coping.



Finished casting with black retentive caps in housing.

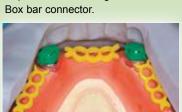

Position the OT Box Classic or Special housings over the rings. Complete the reinforcement using the connectors and join the pieces together with selfpolymerising resin.



Complete prosthesis reinforcement. with cast



Finished wax-up with sprue; ready to be invested.



For additional reinforcement...with the silicon mask in position, insert a wax pin to support each tooth before casting.



### **OT REVERSE 3**




# CLINIC





#### PROSTHESIS WITH REINFORCEMENT IN CAST METAL

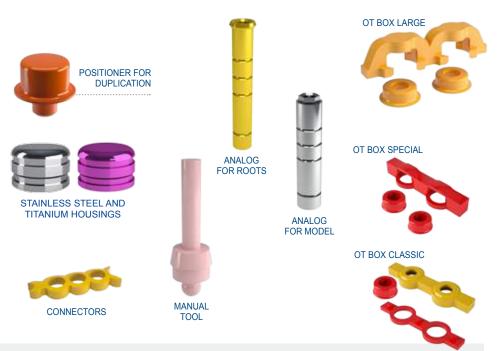

OT REVERSE 3 is a root supported direct pivot attachment system which provides retention and stabillity for full dentures. The "split" male portion of the attachment is manufactured from titanium that is embedded into a soft nylon material. The female pivots have a unique shape that is designed to fit most remaining root structures. OT REVERSE 3 is successful even with minimal bone support of the remaining dentition. The system is cost effective with simple laboratory and chairside procedures.



N. 2 Protective Disks



#### ROOT PREPARATION AND IMPRESSION





Use the diamond sizing to prepare the root for the attachment. Using the hand tool, insert the plastic pivot and apply cement.



Pivots cemented into the relationship in the line into the relationship in the relatio into the pivot and take the impression. For best results, use a stiff bodied impression



The laboratory will place the analog and pour the stone model.



The stone model with the OT REVERSE 3 analog in position.

#### CHAIRSIDE PROCEDURES



If you are using the plastic retentive male, remove the

Caution: If the prosthesis is inserted incorrectly, it could bend and it will not fit into the female housing.



attachment with self-curing resin. It is important to always use the protective disk around the perimeter of the attachment.



When OT Box Large is used, enlarge the space using a carbide bur to reduce interference with the male.



Fill the spaces with selfcuring resin. Insert the prosthesis into the patient's mouth and have them bite down until the resin has cured.



Remove the prostheses and trim the excess resin.

#### FABRICATION OF FRAME FOR DIRECT ROOTS OR IMPLANTS



OT BOX CLASSIC Glue the two OT Box sections together.



OT BOX CLASSIC Separate the two housings and trim any excess material. Use only the part that is needed.



OT BOX SPECIAL Separate the two housings and use only the part that is needed



OT BOX LARGE Separate the two housings and use only the part that is needed.



OT REVERSE 3 Stone model with analogs denture setup and silicon quide.



Insert positioners in the analogs. Apply wax on the gingival crest. Make holes in the wax in contact with the stone. Be sure to use stone separator.



Position the sectioned OT Box housing of choice. Complete the reinforcement using the castable connectors.



Join all of the components with self-curing resin. With the silicon mask in place, in-sert a wax pin for each tooth for additional support.



Remove the OT Box frame from model. Fill in any voids with wax.




Sandblasted Cast Reinforcement



White or pink opaque can be used to block out the metal frame.



The finished prosthesis. Attachments are inserted into the cast housings.

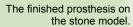




### **RECONSTRUCTIVE SPHERES - OT EQUATOR**

### Titanium + TiN coating



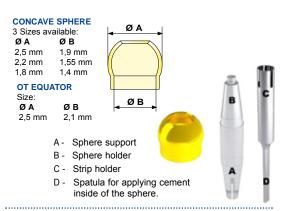

For existing cases with worn spherical attachments which no longer provide adequate retention, the DR8 UNDERSIZED CAP can be used in the early stages of wear of the male component. This elastic cap offers an inner dimension of 1.7 mm and 2.2 mm which is smaller than Rhein83 normal and micro size caps and can be used with standard Rhein83 stainless steel housings.

When ball attachments show excessive wear, the CONCAVE RECONSTRUCTIVE SPHERES are recommended as the best long term restorative option. The CONCAVE RECONSTRUCTIVE SPHERES restore the worn male to it's original size of 1.8 mm, 2.2 mm or 2.5 mm diameter. CONCAVE RECONSTRUCTIVE SPHERES are manufactured with a Titanium Nitrite coating and are rated over 1600 Vickers hard.

The chairside procedure for using the reconstructive spheres is fast, easy and provides an economical alternative to replacing the old restoration.



DR8 Undersized Caps are available in 3 levels of retention for normal and 2 levels of retention for the micro size.

### CLINIC

Dental attachments, like most other mechanisms, are subject to wear out. Rhein83 produces spheres for restoring worn ball attachments which restore and stabalize the prosthesis in a single appointment. Reconstructive spheres are bonded over the worn ball restoring the attachment to it's original size.

# CONCAVE RECONSTRUCTIVE SPHERE RESTORING A WORN OUT SPHERE



Insert the concave sphere into side A of the plastic tool. Fit over the worn out sphere in the mouth.



If the concave sphere does not fit passively, use a cylindrical bur (diamond or carbide) to slightly reduce the diameter. Check the fit again and repeat as needed.



Check the position of the concave sphere on the worn out sphere and finish by cleaning the two parts.



Additional surface can be removed by using side C of the tool. Insert a diamond strip into the notches, place the tool over the sphere and turn the manually.



Place a small amount of two-part self curing "metal to metal" resin inside the sphere.



Place the concave sphere over the worn sphere and wait for the resin to cure.



Once the resin has cured, remove any excess material.



The completed repair. The cap can be repositioned if necessary.



### RECONSTRUCTIVE SPHERES

Titanium + TiN coating

Rhein83 offers two types of reconstructive spheres; A solid sphere and a concave sphere. Both types are titanium nitrate coated with a Vickers hardness rated over 1600. The Concave Reconstructive Spheres are available in 1.8 mm, 2.2 mm and 2.5 mm ball diameter. The Solid Reconstructive Spheres are only available with a 1.8 mm ball diameter. The Concave Sphere is used for restoring worn ball attachments and the Solid Sphere is used for restoring ERA <sup>®</sup> and CEKA <sup>®</sup> type attachments.

#### SOLID RECONSTRUCTIVE SPHERE RESTORING A WORN OUT RING ATTACHMENT



The worn-out female ring attachment.



Apply a small amount of two-part self curing "metal to metal" resin on the bottom of the sphere. Insert the sphere into the attachment using the tool. Wait for the resin to cure.



The female attachment was converted into a male OT Cap Micro directly in the patient's mouth.


#### **MULTIUSE**

SOLID "RECONSTRUCTIVE" TITANIUM + **TIN COATING** rated over 1600 Vickers

TO REBUILD ANY "RING" TYPE **ATTACHMENT** SUCH AS: ERA ® AND CEKA®





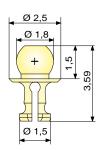





#### SOLID RECONSTRUCTIVE SPHERE RESTORING A WORN OUT OVERDENTURE BAR



Create a hole in the wall of the bar using a 1.6 mm ball drill.



Apply a two part composite to the shank of the sphere. Using the tool, insert the sphere into the hole. Wait for the composite to cure.



The sphere firmly cemented in place. The OT Strategy Cap can now be used in the prosthesis resulting in stability and retention.

#### **OT CEM** is a self and photo curing cement. It is designed for permanent metal to metal bonding in the use of attachments in prosthetic implant solutions. Recommended for the following products:



COPING COVER THREADED SPHERICAL ATTACHMENTS WITH THREADED SLEEVES



#### SOLID RECONSTRUCTIVE SPHERE RECOVERY OF TITANIUM ABUTMENTS



case with unknown titanium abutments. Worn out openings are present on top of the fixtures.



Reconstructive Spheres are placed into the openings. A two-part self curing "metal to metal" resin is applied.



Retentive caps positioned into the existing denture. The denture is now stable and secure.

The SOLID RECONSTRUCTIVE SPHERES can be bonded to the inside of hollow attachments or those with a female ring such as ERA® and


Reconstructive Spheres can be used to repair various attachments available on the market. These attachments can be found in many types of prosthesis including overdentures, implants, roots and frameworks. If worn out or broken, they cannot be repaired easily.

The SOLID RECONSTRUCTIVE SPHERES offer a fast and easy cost effective alternative, transforming a female ring attachment into a male Micro OT CAP attachment. This repair can be completed chairside in a single appointment.



### **BAR AND CAST OVERSTRUCTURE**

on the master model without duplication

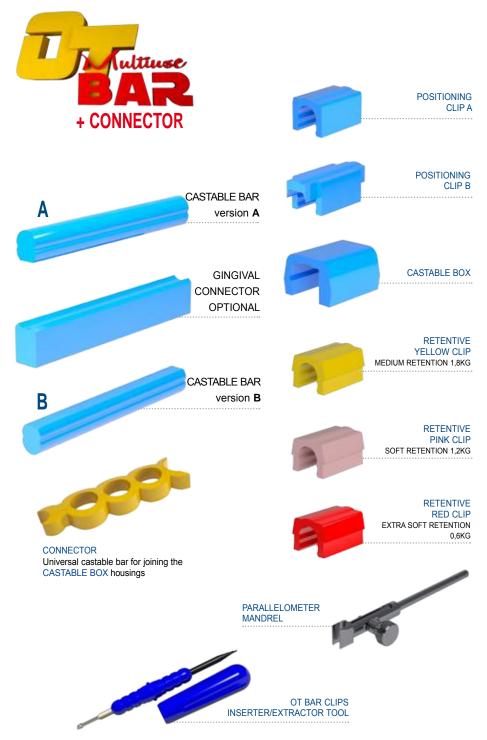

The OT BAR MULTIUSE is designed with a 4 point retentive system. This unique system provides superior retention and can be utilized for both rigid and resilient functionality. With it's innovative two-sided design (Side A is rounded and Side B side is flat), depending on the indication, either side can be used. If a resilient solution is required the bar is positioned with the flat side facing up or if a rigid solution is required then the bar is positioned with the round side facing up. OT BAR MULTIUSE can also be used as a connecting bar between canines in the anterior region.

OT BAR MULTIUSE and the cast housing are fabricated directly on the master model saving time by eliminating the need for duplication.



### LABORATORY










Retention





The rigid bar is used as a "connection" between two stable teeth where a "back and forth" motion is required. The bar can also be used in scenarios involving multiple abutments where the prosthesis is supported by a thin layer of soft tissue.

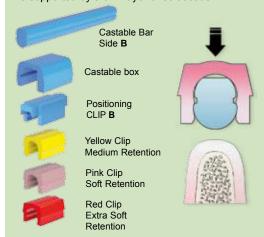

#### SIDE B

The resilient bar is most often used in scenarios involving multuple abutments where the prosthesis is supported by a "normal" layer of soft tissue.

# **CASTABLE BAR IN TWO VERSIONS**

### **RESILIENT - RIGID**















#### FABRICATION OF THE SUPERSTRUCTURE ON THE MASTER MODEL WITHOUT DUPLICATION SIDE A - RIGID



Mount the bar using Side A of the mandrel. Using resin or wax, complete the model.



The finished casting. Be careful not to wear out the retentive surfaces when



Block out any undercuts using wax and place Positioning Clips A on the bar.



To isolate, apply a small piece of tape (ex: teflon, Scotch) on the Positioning Clips A and on the cast bar. Insert the castable box housings.



To prevent resin from adhering to the bar, place a small piece of adhesive tape (ex: teflon, Scotch tape) over the bar. Use self-curing resin to connect the castable boxes.



Complete the model using wax and add castable connectors for extra reinforcement of acrylic. Sprue the model and cast.



The completed casting with retentive clips snapped in place.



The finished denture with cast reinforcment retentive clips in place.

#### SIDE B - RESILIENT



Mount the bar using Side B of the mandrel. Using resin or wax, complete the model.



The completed casting. Use caution when polishing the surface. Be sure not to wear out the retentive undercuts.



Use wax to remove all undercuts. Apply a thin layer of wax on the top of the bar to create a cushion. Insert Positioning Clips B.



To isolate, apply a small piece of tape (ex: teflon, Scotch) on the Positioning Clips B and on the cast bar. Insert the castable box housings.



To prevent resin from adhering to the bar, place a small piece of adhesive tape (ex: teflon, Scotch tape) over the bar. Use self-curing resin to connect the castable



Complete the model using wax and add castable connectors for extra reinforcement of acrylic. Sprue the model and cast.




The completed casting with retentive clips snapped in place.



The finished denture with cast reinforcment and retentive clips in place.




### **EXTRACORONAL CASTABLE ATTACHMENTS**

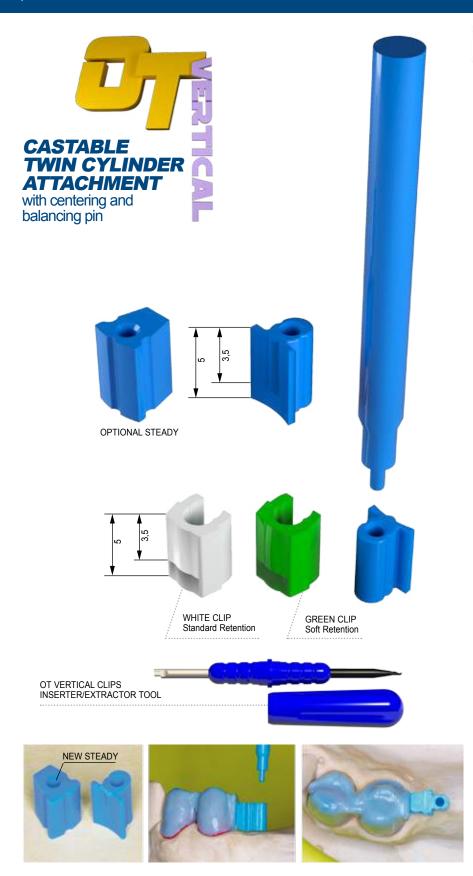

The cast metal guide pin is necessary to center, connect and balance the prosthesis during the final insertion. When milling or "cross arch" stabilization are not possible, the guide pin along with the NEW STEADY will provide lateral stability to the prosthesis. This ensures a longer life for the retentive clips. The vertical height of the attachment can be adjusted by reducing both male and female parts from the original length of 5 mm down to 3.5 mm according to the pre-marked notches. Reducing the vertical height creates no difference in functionality. Removal and replacement of clips can be easily performed by the Dentist chairside.

## **LABORATORY**



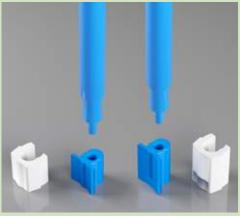



The ot vertical attachments and retentive clips can be lowered in height, in order not to risk inconveniences, it is advisable not to shorten them beyond the indicative sign \* of 3.5 mm.





### REPLACEMENT OF RETENTIVE CAPS

#### TECHNIQUE: SINGLE CASTABLE MALE ASSEMBLY



Parallelometer Key: insert the plug into the hole of the attachment, rotate a few times to obtain the correct seal and be able to slide it out comfortably after fixing the attachment in the wax.



Completed the assembly and the wax modeling, close the hole with wax and create a lingual drilling, proceed with the sprue of the product.



Crowns with cast connection, finished and polished milling.



Retentive clip inserted on the cylindrical male. Waxed model (be careful not to dirty the coping with wax before duplication).



Duplicate model in coating with reproduction of the clip format.



Melted and sanded framework.



Finished work, green retentive clip inserted into the framework.



Work finished.

#### TECHNIQUE: CALSTABLE MALE + STEADY ASSEMBLY



Parallelometer key: insert the plug into the hole of the adapted Steady to follow the mucous profile, rotate to remove the key after fixing it.



Place the Ot Vertical male gluing it lingually to the Steady and complete the waxing.



Crowns with molten attachment, the Steady and Vertical have been adapted to the mucosal profile.



Retentive clip inserted on the cylindrical male. Model discharged with wax (be careful not to dirty the coping with wax before duplication).



Duplicate model in coating with reproduction of the Steady format and clip.



Framework molded, cover with the wax the Steady portion to have an insertion guide and stability when melted.



Framework placed on the Steady



Framework on the model. Stability is guaranteed even without milling, thanks to the Steady

### **COMBINED RETENTION ATTACHMENT**

### For Multi-Functional Prosthetics





**ATTACHMENT** 

#### RETENTIVE CAPS OT CAP Micro



Clear • Standard



Pink • Soft



Yellow • Extra Soft



Green • Elastic



Black • Processing





Clear • Standard





Pink • Soft



Yellow • Extra Soft



Black • Processing

### **TOOLS**







The OT UNILATERAL castable attachment from Rhein83 is specifically intended for unilateral, bilateral or implant bar applications without additional support from milled bracing arms.

OT UNILATERAL's exclusive design features a two-in-one combination of 1.8 mm horizontal and vertical spheres utilizing OT CAP and OT STRATEGY micro size female caps. The male section of the attachment is engineered with a vertical strut which extends through the base of the attachment providing exceptional lateral stability and distal support to the prosthesis.



The Uni-Box female component is a one piece castable housing that covers the entire male section, adding superior strength to the acrylic.

## **LABORATORY**



#### **EXCLUSIVE FEATURE**

2-IN-1 DESIGN - A COMBINATION OF HORIZONTAL AND VERTICAL MICRO SPHERES ARE USED WITH THE OT CAP AND OT STRATEGY ATTACHMENT SYSTEMS





#### **MULTIPLE BENEFITS**

BECAUSE OF IT'S UNIQUE DESIGN, OT UNILATERAL PROVIDES:

- \* LATERAL STABILITY
- \* NO MILLING REQUIRED
- \* SUPERIOR RETENTION
- \* CONTROLLED RESILIENCY
- \* OVERALL FUNCTIONALITY
- \* ECONOMICAL SOLUTIONS



### Superstructure Set Up Technique

#### UNILATERAL SADDLE: ATTACHMENT AND OVERSTRUCTURE UNIQUE PHASE SETTING UP



Positioning of the OT UNILATERAL bar using the OT CAP paralleling mandrel by starting with the analysis of the masticatory plan. Proceed by connecting the bar to the last modeled wax crown.



Place the positioning ring over the OT CAP micro sphere. Place the castable OT BOX component in position, the positioning ring will assure the proper position.



Join the Uni-Box component to the connector by using a pattern resin in order to reinforce the structure. Be careful not to have any material inside the Uni-Box component.



Remove the positioning ring by the OT CAP sphere and proceed with the sprue procedure.



Unique fusion is one of the best features of the UNILATERAL attachment.



Fused UNILATERAL and Uni-Box. Sandblast the casting by keeping attention not to "over-sandblast" the spheres. Insert the black laboratory caps and proceed by polishing the sphere.



In order to provide the optimal stability. wax-up carefully the saddle in order to embrace the ridge as much as possible.



Completed procedure: proper retentive caps (adeguate degrees of elasticity) are placed inside the fused Uni-Box component

#### BILATERAL STRUCTURE: RESILIENT FUNCTIONALITY AND FREE MILLING PROCEDURE



Place the positioning ring over the OT CAP micro sphere. Place the castable OT BOX component in position, the positioning ring will assure the proper position.



Finished work: Ot cap and Ot Strategy caps, with the proper retention features, are inserted inside the Ot-Box component.

#### **IMPLANT SUPPORTED BAR:** DISTAL EXTENSIONS AND COMBINED FUNCTIONALITY



Once the components to build the bar are inserted, place the OT UNILATERAL bar by using the OT CAP mandrel and by analyzing the masticatory plan. Connect it then distally to the modeled bar.




Cast bar thank to the combined functionality of the OT UNILATERAL. The prosthesis will count on a improved stability without any additional stress over the implants.



### **LOCKING PIN - TITANIUM**



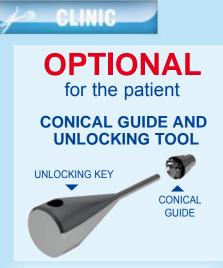




UNLOCK THE PROSTHESIS BY INSERTING THE UNLOCKING TOOL INTO THE CONICAL GUIDE.


### ADJUSTABLE TITANIUM LOCKING PIN


SPACER RING SYSTEM TO POSITION THE KEY TO THE DESIRED SHAPE



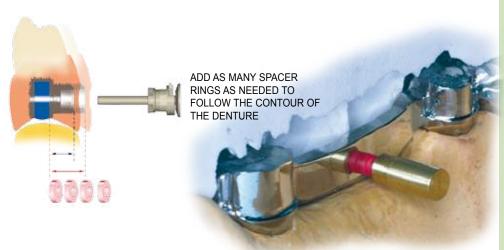

## LABORATORY

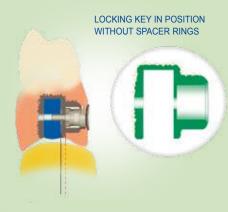

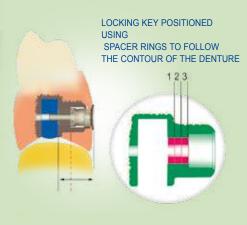

### **LOCKING PIN - TITANIUM**













Model the bar in resin and drill a 0.8 mm hole in the most ideal position.



Insert the ceramic pin through the



The finished and polished bar.



Insert the housing shaper into the hole and lock it in place using resin. Be sure not go past the "STOP" when appling resin



Using resin, complete the model of the superstructure up to the "STOP". Remove the housing shaper and



Pull out the brass positioner and cast.



Insert the pre-fabricated housing and bond.



Insert the positioner again. Proceed with wax and cure the



Insert the locking key into the prefabricted housing guide. The "keyring" mechanism is now locked.



Bend the locking key and brake



Apply the self-hardener composite material to stop the locking key and insert the locking pin in the hole.



Locking Pin locked in position. Finish and polish.





Finished prosthesis. Determine whether or not to use the EXTRACTOR KEY



### **OVERDENTURE ATTACHMENTS - SPHERO FLEX - SPHERO BLOCK**

Rotating & Stationary Ball Abutments For Divergence Correction



The Sphero Flex implant overdenture attachment is compatible with all implant systems currently on the market. Featuring a rotating ball with a diameter of 2.5 mm that is flexible to 7.5° in all directions. When used with a 14° directional ring, Sphero Flex corrects divergence up to 43° between two implants. Sphero Flex creates a passive path of insertion which reduces trauma to the implant.

Sphero Block is a "one-piece" milled stationary ball implant attachment. It is available in 2.5 mm and 1.8 mm diameters. Sphero Block provides exceptional stability and corrects divergence up to 28° between 2 implants Sphero Block implant attachments are compatible with all implant systems currently on the market.

Sphero Flex and Sphero Block are manufactured with cuff heights ranging from 1 mm to 7 NOTE: The Sphero Flex and Sphero Block attachments are available for all platform diameters.

## GLINIC CLINIC







#### DIRECTIONAL RINGS CORRECT PLACEMENT

Before placing the impression abutment on the implant it is suggested to put a gray directional ring (for parallel systems) or a ring for angled implants if not parallel. This will keep the impression coping "on level" during the impression. The directional rings have only one direction of insertion.



**LABORATORY** 





3 EASY STEPS

1. Place directional rings (green and red are shown here) over the spheres establishing a level plane.

- 2. OT BOX positioners are placed over spheres to support box housing during framework fabrication.
- 3. After gluing the 2 OT BOX parts, cut and use the necessary pieces for the housing.

### **OVERDENTURE ATTACHMENTS - SPHERO FLEX - SPHERO BLOCK**

Rotating & Stationary Ball Abutments For Divergence Correction

#### CHAIRSIDE PROCEDURE FOR POSITIONING THE CAPS



Screw the attachment into the implant. For best results, unscrew and screw the attachment 3/a times and then tight firmly.



Select the appropriate directional rings and place them over the spheres. Be sure that the ring is aligned with the hex and seated properly on the platform.



Once the directional rings have been positioned, it is advisable to remove the rententive caps and place a protective disk over the spheres. Replace the retentive caps in original position when finished.



Try the prosthesis in the mouth. Check to see if there is enough space for the retentive caps. Fill the holes with self-curing resin and position the prosthesis over the caps and spheres in the patient's mouth.



Once the resin has hardened, remove the prosthesis. Remove the protective disk along with any excess resin.



Finished prosthesis

#### TAKING IMPRESSION TRANSFER



Place the directional ring over the sphere with the flat side facing down. Place the impression coping over the sphere.



Rotate the directional rings to achieve a common axis parallel to the occlusal plane and take the impression.



Remove impression. Directional rings must be removed from the impression and spheres.



Place the analogs into the impression copings and send to the laboratory for model fabrication.

### OT BOX CLASSIC NORMAL - CAST REINFORCED ACRYLIC PROSTHESIS USING DIRECTIONAL RINGS



Place directional rings over the spheres. OT BOX is placed over the directional rings, ensuring that the horizontal plane is level. Connect with resin.



OT BOX The constructed substructure with reinforced wax pins. Sprued and ready for casting.



cast substructure on the model. The metal reinforcement pins for each tooth are positioned according to the silicone mask.



Finished prosthesis with caps inserted in the cast OT BOX housings.



Screw the abutment into the analog. Be sure to use the abutment with the proper cuff height.



Directional rings are placed over the abutments and must be fully seated on the platform. Rotate rings until they are parallel in the same horizontal plane.



The nylon caps are inserted into the stainless steel housings and placed on top of the directional rings. Verify that the caps are still in the same horizontal plane.



The finished prosthesis with stainless steel housings and retentive caps in final position.




### **IMPLANT OVERDENTURE ATTACHMENTS**

**Components and Accessories** 

#### DIRECTIONAL RINGS

for angle correction

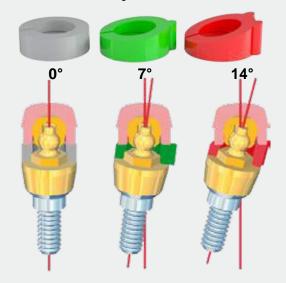

#### SPHERO FLEX - SPHERO BLOCK

In order to achieve a passive fit for the final prosthesis using the SPHERO FLEX and SPHERO BLOCK attachment systems, it is necessary to use DIRECTIONAL RINGS. When not used, there is a high possibility that the attachments will not seat properly into the prosthesis due to incorrect positioning of the caps. This mis-alignment will result in premature wear of the caps causing additional trauma to the implant. SELECTION OF DIRECTIONAL RINGS: The position and angulation of the implant will determine which directional ring will be used. For parallel implants, a 0° DIRECTIONAL RING can be used. For implants that have greater divergence, a 7° or 14° ring can be used.

Place the DIRECTIONAL RING onto the hex of the attachment with the flat side down. Be sure that the ring is fully seated. Next, place the retentive cap onto the sphere and rotate the DIRECTIONAL RING until the cap is parallel with the other caps and are in the same horizontal plane. This ensures that the retentive caps are correctly alligned inside of the final prosthesis.



## INSTRUCTIONS FOR USE OF ABUTMENT DRIVER / WRENCH

Abutment Driver has a sliding mechanism that locks it onto the ball abutment. This needs to be fully engaged to properly tight the abutment without damaging the abutment. To dis-engage driver once the abutment is tightened in the mouth push down on the silver portion to loosen the driver from the abutment (Please screw and unscrew the abutment 3/4 times in order to achieve a fine adaption of the two threads). Then tight the abutment with a torque controller or the manual torque wrench.









Clamping mechanism

Incompletely seated driver

Driver fully seated

# UNIVERSAL "ANTI-UNSCREWING" SYSTEM WITH ELASTIC INSERT Recommended for Sphero Flex, Sphero Block and

Recommended for Sphero Flex, Sphero Block and OT Equator attachments with a cuff height over 5 mm. This system can also be used for single screws. (Core Vent, Branemark, Pitt Easy, Bona Fit)





#### **ELASTIC INSERT**

This component is manufacutred from bio-compatible materials with an "elastic" memory. While screwing in the attachment, the insert is compressed. When the threaded attachment is fully seated, the elastic insert will expand and return to it's original form, which prevents rotation and unscrewing of the device. The insert is applied at the manufacturing facility UPON REQUEST. It can be applied to any screw with a diameter greater than 1.8 mm.



# MINI PARALLELOMETER WITH MODEL HOLDER BASE

### LABORATOR

### **FEATURES:**

- EASY TO USE
- COMPACT
- PRECISE
- ECONOMICAL

Height of the mini-parallelometer: 17cm



divergence indicator



The MINI-PARALLELOMETER allows accurate positioning of attachments without the need for an expensive milling machine. The MINI-PARALLELOMETER is a useful and economical device for the laboratory technician that can be used in day-to-day operations or in a training environment.

#### INSTRUCTIONS FOR USE

Place the stone model on the swivel base. Rotate the base until the ideal model position is found. Insert the mandrel into the notch on the horizontal extension arm and lock it into place by tightening the screw. Adjust the height by moving the horizontal arm up and down. Once the correct height has been found, lock the arm into position by tightening the rear locking screw.



To determine the tissue height above the implant and eliminate mistakes when choosing the correct attachment, the Cuff Height Measuring tool is reccomended. The Cuff Height Tool is compatible with all implants that have an internal or external hex connection.

#### **CUFF HEIGHT MEASURING TOOL FOR IMPLANTS**

- 1. Rotate upwards the gold colored plate until the tool is completely open.
- 2. Insert the tool into the implant. Be sure that it is fully seated on the top of the implant.
- 3. Firmly hold the tool and rotate the gold plate clockwise until it contacts the ridge.
- 4. Remove the tool and read the color coded rings indicated on the pin to determine the cuff height.

When a colored ring is completely covered, and only the silver band between colors is visible, it is recommended to utilize the next (higher) color.

#### IMPORTANT:


Before ordering an attachment, it is necessary to specify: Implant manufacturer, implant brand, diameter, internal or external hex connection and cuff height. The cuff height is determined by taking the corresponding color from the cuff height measuring tool. For implants with an internal hex connection the cuff height will range from .5 mm to 7 mm and for implants with an external hex connection, the cuff height will range from 1 mm to 7 mm.

#### **Cuff Height Measuring Tool** With Stationary Pin

The cuff height measuring tool with stationary pin provides the same functionality as the tool with a threaded pin, however it is used in cases where there is limited space between two implants.

#### **Cuff Height Measuring Tool** With Threaded Pin And Ball Indicator

With easy to read color-coded millimeter measurements. Dentists and dental laboratories can accurately measure tissue height between 0.5 mm and 7 mm. The ball indicator outlines where the male component of the attachment will seat above the tissue.

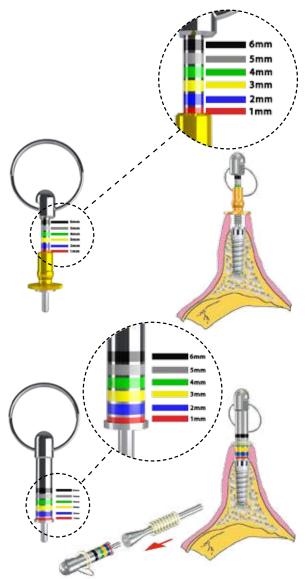

## **BROKEN SCREW EXTRACTOR KIT FOR IMPLANTS**

FOR REMOVAL OF BROKEN IMPLANT SCREWS



inserted the CLAW REAMER BUR (C) in the MANUAL DEVICE (B) for the manual removal of the broken screw





## COMPONENTS AND ACCESSORIES:

- A CENTERING DEVICE
- **B MANUAL DEVICE**
- C CLAW REAMER BUR
- **D REVERSE CUTTING BUR**

### **BROKEN SCREW EXTRACTOR KIT FOR IMPLANTS**

FOR REMOVAL OF BROKEN IMPLANT SCREWS

A broken screw inside an implant is a serious, even if not so frequent, problem. With the BROKEN SCREW EXTRACTOR KIT, you can remove the broken screw fragment from the implant if it has not been cemented or if the implant internal thread has not been damaged in a previous removal attempt. In 90% of the cases the broken screw can be easily unscrew but, the operation must be carried out with great skill, patience and attention. The time necessary for the removal may depend on a number of factors, including the location of the implant which may facilitate or complicate the operation. Once the screw has been removed, the leftovers must be removed from the implant with air, water, and suction.


#### **WARNING:**

During the use it is mandatory to cool down the CENTERING DEVICE (A), the CLAW REAMER BUR (C) and the REVERSE CUTTING BUR (D) with a lot of water in order to not overheat the implant; consequently, the bone will be protected from any risk of overheating and necrosis. The effectiveness of the CLAW REAMER BUR and of the REVERSE CUTTING BUR is optimal for three extractions of broken screws. The REVERSE CUTTING BUR is extremely hard but brittle to bending; in order to avoid its breaking it is absolutely necessary to that the CENTERING DEVICE does not move during the entire operation. For some types of connection the BROKEN SCREW EXTRACTOR KIT is available in stock; for other types of connection it is necessary to start production and the production time increases to a maximum of 10 working days.

#### **INSTRUCTIONS FOR USE:**

Insert the CLAW REAMER BUR into the CENTERING DEVICE by inserting them into the implant and with a constant pressure on the broken screw through the MANUAL DEVICE (B), unscrew it and remove it from the implant (Fig.1 and 2). In case the screw remains locked inside the implant, insert the CLAW REAMER BUR into the 20:1 implantology contra-angle handpiece. Set the implant insertion program in anti-clockwise rotation and always through the MANUAL DEVICE (B) insist on the broken screw to flatten it and prepare it for the REVERSE CUTTING BUR that will destroy it (Fig.3). Insert the REVERSE CUTTING BUR in the 20:1 implantology contra-angle. Set the program in anticlockwise rotation with speed between 500 and 600 rpm. Refrigerate with plenty of water during this operation. Insert the REVERSE CUTTING BUR into the CENTERING DEVICE, start the rotation, keep it pressed for no more than 3 seconds on the broken screw and release. This alternating movement facilitates the entry of water for the refrigeration of the implant and of the REVERSE CUTTING BUR. It is absolutely necessary that the CENTERING DEVICE does not move during the whole operation (Fig.4). If the CENTERING DEVICE moves, the REVERSE CUTTING BUR will break. Once all laser engravings of the REVERSE CUTTING BUR disappear in the CENTERING DEVICE the operation is complete and the screw is completely destroyed.

Once the screw has been destroyes, the leftovers must be removed from the implant with air, water, and suction.

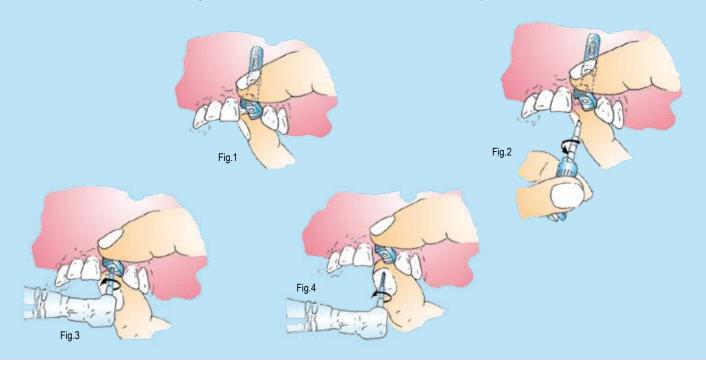

### **POLISHING BURS KIT BY CARLO BORROMEO**



#### **METAL POLISHING KIT**

From the experience of Carlo Borromeo, dental technician with thirty years of experience in metallurgy and casting, precursor of the use of digital techniques for the construction of metal structures, has come to light a set dedicated to the metal polishing. The aim was to identify a simple protocol that could be applied by everyone.

Well polished metals surfaces avoid or delay the plaque adhesion. Clean surfaces are very important to reduce the risk of perimplantitis. Also the retention systems (attachments) should be carefully polished so to allow them to work correctly and for much longer reducing the possibility of abrasion and therefore of malfunction.

The internal walls of the counter-bar if well polished prevent the formation of plaque, thus increasing the durability of the attachments and the functionality of the denture.









Polishing milled parts



Polishing extended parts



Finishing





Polishing small parts around connections





Rubbers with fine diamond grains







### **INSTRUCTIONS AND TECHNICAL ADVICE**



#### REPLACEMENT OF CAPS

Rhein83 recommends that caps should be replaced every 12 months. The longevity of the caps is affected by many variables including: original case design, patient hygiene and general maintenance of the prosthesis.



#### HOW TO REPLACE THE CAPS

In a prosthesis with metal housings, the cap can be removed by using the extractor tool for caps; otherwise use a spherical bur at low RPM without damaging the housing.



In a prosthesis where the cap is incorporated directly into the resin, it can be removed by hand with a pointed instrument (such as a spatula) or the Rhein83 cap extractor tool. If a bur is used, be careful to remove only the retentive cap and to not modify the form that remains in the resin. If the resin site is damaged during the removal of the cap, repair the area with self curing resin before inserting the new cap. The cap insertion tool is used for this procedure.



#### **GREEN ELASTIC CAPS**

These caps are highly elastic and have a medium level of retention. In cases where metal housings are used, it is recommended to apply a drop of adhesive (cyanoacrylic) on the inside of the housing before inserting the cap.



#### **TITAN CAPS**

These caps were designed to be used on the OT CAP TECNO as well as the Normal and Micro attachments with machined titanium spheres.



#### **CAP INSERTION TOOL**

When using high retention caps, it is recommended to insert them directly in the clinic into the housing using the cap insertion tool.

OT CAP Normal / Micro OT Reverse.



#### PROSTHESES WITH MULTIPLE ATTACHMENTS

In order to balance the retentive levels of a prosthesis with multiple attachments, it is possible to use caps with different levels of retention.

REAMERS AND CAP TESTERS: if the retention of the caps is too hight, insert the reamer into the caps and rotate it in a clockwise direction, after only a few rotations it will wear down the perimeter which will reduce the retention. Try the prosthesis in the mouth, if it is still too retentive, repeat the operation with the reamer, In order to avoid trying the prosthesis in the mouth too many times, one can use the spherical tester, in order to evaluate the holding strength.



#### HOW THE RETENTIVE CAP FUNCTIONS

The Rhein83 caps are manufactured with a high elasticity which creates both mechanical and frictional retention resulting in a larger contact zone between the cap and the lower portion of the sphere. A small space between the metal housing and the cap allows the cap to expand as it passes over the equator of the sphere. Once completely engaged, the cap returns to its original form.



POLISHING OF THE "CAST" ATTACHMENTS: It is recommended that only glass beads or a soft cloth wheel are used to polish attachments. In order to avoid damage to the sphere duing these procedures, it is a good practice to cover the spheres with a retentive cap. The retentive caps can be reused again for this procedure.



### **DEMONSTRATION MODELS**

TRADI	TIONAL
PROS1	THESES





06P Model with upper prosthesis with OT Cap Normal / Micro size attachments:

1 OT CAP NORMAL 1 OT CAP MICRO

1 Frame with OT BOX mono housings

5 Acrylic teeth



07P Model with lower prosthesis with OT Strategy

1 OT STRATEGY

1 OT STRATEGY + STEADY

1 Frame with caps and duplicated housings

5 Acrylic teeth



PROSTHESIS ON NATURAL TEETH

Model with lower "Overdenture" prosthesis:

1 PIVOT FLEX titanium post1 Cast post with OT CAP sphere1 Complete denture with 14 teeth

1 Cast OT BOX reinforcement incorporated in the denture



04P/A Same model as 04P. Denture with pre-fabricated STAINLESS STEEL HOUSINGS for retentive caps

OTAINEEGO OTELE HODOINGO IOI Telentive caps

09P

04P

MODEL WITH LOWER PROSTHESIS WITH OT VERTICAL

1 OT VERTICAL

1 OT VERTICAL + STEADY

1 Frame with clips and duplicated housing

6 Acrylic teeth



IMPLANT MODELS

Model with lower prosthesis with SPHERO FLEX abutments:

2 Implant analogs 1 SPHERO FLEX

1 SPHERO BLOCK

1 Complete denture with 14 teeth

1 Cast OT BOX reinforcement incorporated in the denture

031/A

08B

031

Same model as 031. Denture with pre-fabricated

STAINLESS STEEL HOUSINGS



PROSTHESIS ON FIXTURES

Model with lower prosthesis with OT Bar Multiuse:

2 Implant analogs


1 Cast bar with copings

1 OT BAR MULTIUSE


1 Cast superstructure with two retentive clips


1 Complete denture with 14 teeth

### **PRODUCT RANGE - SIZES AND DIMENSIONS**

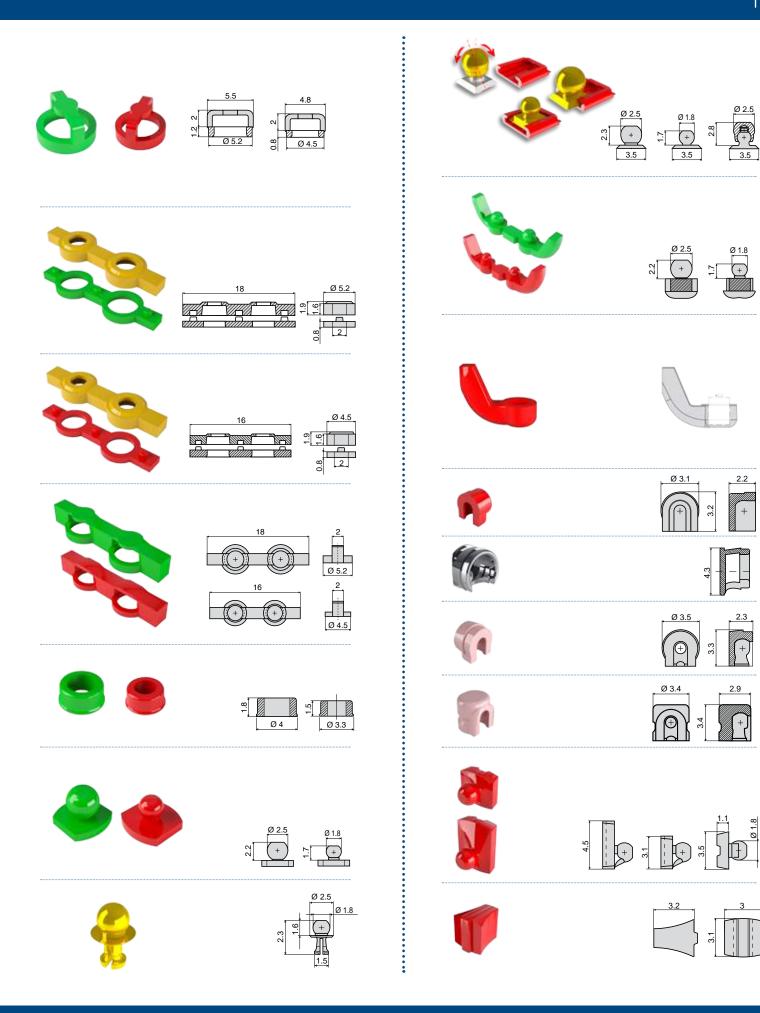

## **PRODUCT RANGE - SIZES AND DIMENSIONS**

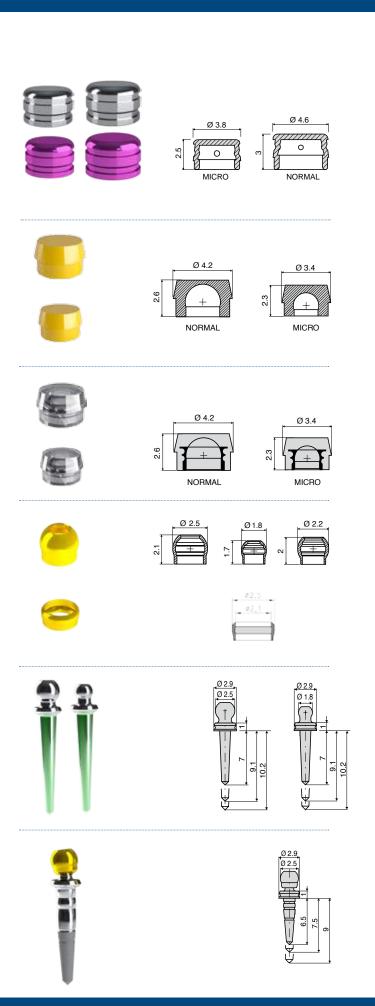


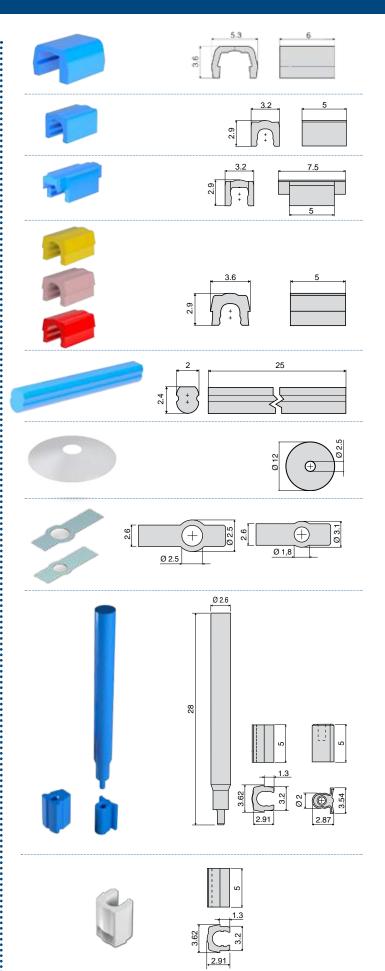

## **PRODUCT RANGE - SIZES AND DIMENSIONS**





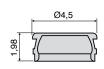























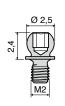









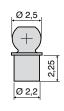



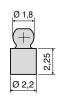




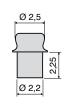








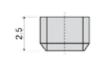









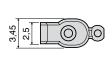





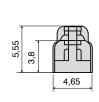

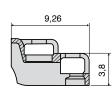

### KITS AND CODES



Ref.: 005SKLUS

INTRODUCTORY KIT FOR LABORATORY contains assorted attachments and tools



#### Ref.: 011SPL

S.P.L. INTRODUCTORY KIT contains assorted titanium BLOCK pivots, castable pivots, retentive caps and two regulating tools


#### ATTACHMENTS OT CAP SYSTEM




Ref.: 038STF

Ref.: 038STN



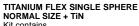

- Kit contains
  2 SINGLE TITANIUM SPHERES
  2 PINK CAPS SOFT RETENTION
  1 SPHERE HOLDER
- CASTABLE SLIDING BASES

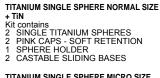

#### TITANIUM SINGLE SPHERE MICRO SIZE




Ref.: 038STM

- + TiN
- Kit contains
  2 SINGLE TITANIUM SPHERES
  2 PINK CAPS SOFT RETENTION
  1 SPHERE HOLDER
- CASTABLE SLIDING BASES



Ref.: 093CTN NORMAL

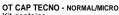

- TRANSPARENT CAPS NORMAL/MICRO

  OTHER TRANSPARENT CAPS NORMAL/MICRO
- PINK CAPS NORMAL/MICRO YELLOW CAPS NORMAL/MICRO



Ref 093CTM MICRO

Ref.:

092CAN

OT CAP NORMAL SIZE Kit contains

- 4 SINGLE SPHERES
- 4 PINK RETENTIVE CAPS
- 4 STAINLESS STEEL HOUSINGS
- (2 for resin 2 for soldering)
- 4 PLASTIC POSITIONING RINGS



Ref.: 092CAM

#### OT CAP MICRO SIZE

- Kit contains SINGLE SPHERES
- 4 PINK RETENTIVE CAPS
- STAINLESS STEEL HOUSINGS (2 for resin - 2 for soldering)
- 4 PLASTIC POSITIONING RINGS



Ref.: 196BCN

#### "ECONOMIC" OT CAP NORMAL SIZE

- Kit contains
- CASTABLE BAR
- CASTABLE BEVELLED BAR
- CLEAR RETENTIVE CAPS
- (Standard retention)



Ref.: 197BCM

#### "ECONOMIC" OT CAP MICRO SIZE

- Kit contains
- CASTABLE BAR
- CASTABLE BEVELLED BAR
- CLEAR RETENTIVE CAPS
- (Standard retention)



Ref.: 099BSN

#### OT CAP & MONO OT BOX FOR FRAME NORMAL SIZE

- Kit contains
- 2 CASTABLE BARS (1 straight 1 bevelled) 4 CLEAR RETENTIVE CAPS
- CASTABLE MONO OT BOX
- PLASTIC POSITIONING RINGS



Ref.: 099BSM

#### OT CAP & MONO OT BOX FOR FRAME MICRO SIZE

- Kit contains
- 2 CASTABLE BARS (1 straight 1 bevelled)
- CLEAR RETENTIVE CAPS
- CASTABLE MONO OT BOX
- 4 PLASTIC POSITIONING RINGS



Ref.: 058BSN

#### OT BOX SPECIAL NORMAL SIZE + CONNECTORS

Kit contains

- 2 OT BOX SPECIAL BARS
- PLASTIC POSITIONERS
- CONNECTORS



058BSM

#### OT BOX SPECIAL MICRO SIZE + CONNECTORS

- Kit contains
- 2 OT BOX SPECIAL BARS
- 4 PLASTIC POSITIONERS
- CONNECTORS



Ref: 153BCN

#### OT BOX CLASSIC NORMAL SIZE + CONNECTORS

- Kit contains
- 2 UPPER BARS
- LOWER BARS PLASTIC POSITIONERS
- CONNECTORS



Ref.: 153BCM

#### OT BOX CLASSIC MICRO SIZE + CONNECTORS

- 2 UPPER BARS
- LOWER BARS
- PLASTIC POSITIONERS
- CONNECTORS



Ref.: 087CRS

#### CONCAVE RECONSTRUCTIVE SPHERE

- Kit contains
- CONCAVE SPHERES IN TITANIUM TIN COATED
- 2 PINK CAPS SOFT RETENTION
- INSERTION TOOL
- GAUGE AND STRIP HOLDER
- Available in 1.8 mm, 2.2 mm, 2.5 mm diameters



Ref.: 087CRE

#### CONCAVE REPAIR OT EQUATOR

- Kit contains
- 2 CONCAVE OT EQUATOR IN TITANIUM TIN COATED
- 2 PINK CAPS SOFT RETENTION
- INSERTION TOOL
- GAUGE AND STRIP HOLDER



Ref.: 089SRS

#### SOLID RECONSTRUCTIVE SPHERE Kit contains

- SOLID SPHERES IN TITANIUM TIN COATED PINK CAPS SOFT RETENTION
- 2 PROTECTIVE DISKS
- KEY
- Available in 1.8 mm diameter



Ref: 064ACN **NORMAL** 

#### ASSORTED RETENTIVE CAPS

- Kit NORMAL Kit MICRO
- 6 CLEAR CAPS STANDARD RETENTION
- 6 PINK CAPS SOFT RETENTION
- 6 YELLOW CAPS EXTRA SOFT RETENTION 6 GREEN CAPS VERY ELASTIC RETENTION



Ref.: 064ACM **MICRO** 

#### ASSORTED RETENTIVE CAPS

- Kit NORMAL Kit MICRO 6 CLEAR CAPS STANDARD RETENTION
- 6 PINK CAPS SOFT RETENTION 6 YELLOW CAPS - EXTRA SOFT RETENTION
- 6 GREEN CAPS VERY ELASTIC RETENTION

#### S.P.L. PIVOTS



**PIVOT FLEX - TITANIUM PIVOT** WITH SWIVEL SPHERE NORMAL SIZE (Ø 2.5 mm)

FOR DIRECT OVERDENTURE

(3 Sizes available) Kit contains

- TITANIUM PIVOT WITH ROTATING SPHERE (adapted for COPING COVER)
- STAINLESS STEEL HOUSING
- FOR RESIN
- 2 PINK CAPS Normal Size Soft retention
- ALUMINIUM DISK
- 3 DIRECTIONAL RINGS



Ref.: 036PTN

**TITANIUM PIVOTS Normal Size** Adapted for COPING COVER Kit contains

5 TITANIUM PIVOTS Sphere 2.5 mm

TITANIUM PIVOTS Micro Size Ref.: Adapted for COPING COVER Kit contains

5 TITANIUM PIVOTS Sphere 1.8 mm

036PTM



Ref.: **Ref. 336PTE7 Ref. 336PTE9** Ref. 336PTE10

OT EQUATOR PIVOTS



Ref : 010PSP

CASTABLE PIVOTS NORMAL SIZE



Ref.:

012PSM CASTABLE PIVOTS MICRO SIZE



Ref.: MOOSER BUR

A01MOG Reamer for post 7 mm, 9 mm, 10 mm



Ref.: MOOSER BUR

A03MOB Reamer for post 12 mm, 14 mm



Ref.: 485IC

OT CAP NORMAL AND MICRO SIZE CAPS INSERTER/EXTRACTOR TOOL



Ref.: 74AC01

PARALLELOMETER MANDREL FOR OT CAP NORMAL SIZE



Ref.: 74AC02

PARALLELOMETER MANDREL FOR OT CAP MICRO SIZE



Ref: 74AC03 PARALLELOMETER MANDREL FOR OT CAP TECNO NORMAL AND MICRO SIZE



Ref.: 080RCN

REAMER TOOL FOR CAPS OT CAP NORMAL SIZE



Ref.: 080RCM

REAMER TOOL FOR CAPS OT CAP MICRO SIZE



Ref: 082ATN

TOOL FOR TESTING CAP **RETENTION OT CAP NORMAL SIZE** 



083ATM

TOOL FOR TESTING CAP RETENTION OT CAP MICRO SIZE



Ref.: 491EC CAPS EXTRACTOR TOOL WITH MULTIUSE HOUSING FOR CLIPS AND CAPS INSERTION

#### **OT STRATEGY ATTACHMENTS**



Ref: 098SSS OT STRATEGY CAPS FOR DUPLICATION TECHNIQUE Kit contains CASTABLE MALES 2 Standard + 2 High

CASTAB:LE STEADY RETENTIVE CAPS



Ref.: 098CAL OT STRATEGY CAPS FOR DUPLICATION TECHNIQUE Kit contains

it contains
CASTABLE MALES (2 Standard + 2 High)
CASTABLE STEADY
STAINLESS STEEL HOUSINGS
POSITIONING RINGS

RETENTIVE CAPS



Ref : 047ACS OT STRATEGY ASSORTMENT CAP KIT FOR DUPLICATION TECNIQUE Kit contains

YELLOW CAPS - EXTRA SOFT RETENTION

4 PINK CAPS - SOFT RETENTION
4 CLEAR CAPS - STANDARD RETENTION



Ref.: 045ACS OT STRATEGY ASSORTMENT CAP KIT FOR STAINLESS STEEL HOUSINGS Kit contains

YELLOW CAPS - EXTRA SOFT RETENTION PINK CAPS - SOFT RETENTION CLEAR CAPS - STANDARD RETENTION

Ref.: 486ICS

OT STRATEGY CAPS INSERTER/EXTRACTOR TOOL



Ref.: 75AC04

PARALLELOMETER MANDREL FOR OT STRATEGY



Ref.:

REAMER TOOL FOR 081RCS OT STRATEGY CAPS

#### OT BAR MULTIUSE ATTACHMENTS



Ref.: 0210BM

OT BAR MULTIUSE Kit contains BARS

POSITIONING CLIPS A POSITIONING CLIPS B

BOXES

RETENTIVE PINK CLIPS RETENTIVE YELLOW CLIPS CONNECTORS **GINGIVAL CONNECTOR** 



Ref.:

OT BAR CLIPS 429IOBM INSERTER/EXTRACTOR TOOL



028OCP

PARALLELOMETER MANDREL FOR OT BAR MULTIUSE

#### OT VERTICAL ATTACHMENTS



Ref.: 0710BV

OT VERTICAL Kit contains

4 CASTABLE MALES CASTABLE STEADY

RETENTIVE WHITE CLIPS RETENTIVE GREEN CLIPS CASTABLE PARALLELOMETER KEYS + PIN



Ref.: 472ICV

OT VERTICAL CLIPS INSERTER/EXTRACTOR TOOL



#### **OT EQUATOR CASTABLE**



OT EQUATOR CASTABLE

Kit Contains:

- 2 CASTABLE MALES 2 TITANIUM HOUSINGS
- **4 RETENTIVE CAPS**

#### OT EQUATOR FOR IMPLANTS

Ref.: 130

Ref

OT EQUATOR FOR IMPLANTS

- Kit Contains:
   1 OT EQUATOR
- TITANIUM HOUSING
- **4 RETENTIVE CAPS**



Ref. 030

1 OT EQUATOR ABUTMENT \*Compatible with ALL implant systems



Ref : 335SBC OT EQUATOR SMARTBOX KIT self-aligning caps housing

Kit Contains:

- 1 SMARTBOX HOUSING WITH BLACK CAP FOR LABORATORY
- 1 PINK PROTECTIVE DISK 4 RETENTIVE CAPS
- - (1 EXTRA-SOFT, 1 SOFT, 1 STANDARD, 1 STRONG)

Ref.: 330SBE Kit Contains:

1 SMARTBOX HOUSING WITH BLACK CAP FOR LABORATORY

#### OT EQUATOR BAR



Ref. 160EQB OT EQUATOR WITH THREADED SLEEVE For Bonding

Kit Contains:

- 2 THREADED OT EQUATOR-1.6 mm thread 2 THREADED SLEEVES-1.6 mm thread
- 2 STAINLESS STEEL HOUSINGS
- WAXING SPACERS
- 8 RETENTIVE CAPS 2 YELLOW EXTRA SOFT

  - 2 YELLOW EXTRA SOFT 2 PINK SOFT 2 CLEAR STANDARD 2 BLACK PROCESSING



Ref. 039SFE2

- 1 THREADED OT EQUATOR
  - 2 mm universal thread

#### OT EQUATOR ELASTIC SEEGER



Passive Bar Connection **ELASTIC SEEGER** 

Ref.: 158ESA Kit Contains:

- t Contains:

  1 CASTABLE CYLINDER HOUSINGS
  FOR SEEGER
  1 SELF-EXTRACTING SEEGER
  1 TITANIUM LOCKING SCREW
  FOR SELF-EXTRACTING SEEGER

#### **TOOLS**



Ref. 74AC01

 PARALLELOMETER MANDREL **NORMAL** 



Ref.: **774CHE** 

- OT EQUATOR SQUARE DRIVER 1.25 mm + HOLDER
- Ref. 760CE
- OT EQUATOR HANDPIECE CONNECTOR 1.25 mm



Ref.: 085818

· STEEL INSERTION TOOL FOR SEEGER



Ref.

CAPS EXTRACTOR TOOL WITH UNIVERSAL INSERTER HOUSING

#### **ACCESSORIES**



▶ 044CAIN • 2 IMPRESSION TRANSFER ( pick up impression )



144MTE • 2 IMPRESSION TRANSFER



Ref.:



144AE

2 STAINLESS STEEL ANALOGS For OT Equator



Ref.: 485IC

CAPS INSERTER/EXTRACTOR TOOL (OT EQUATOR-NORMO-MICRO)

#### SPARE PARTS



Ref.:

- 192ECE
- Kit Contains:

   1 STAINLESS STEEL HOUSING

   1 BLACK CAP PROCESSING

   4 RETENTIVE CAPS:

  1 YELLOW EXTRA SOFT 1 PINK SOFT

  1 CLEAR STANDARD 1 VIOLET RIGID

  1 BLACK PROCESSING 1 PROTECTIVE DISK



#### OT CAP / OT EQUATOR IMPRESSION COPINGS



Ref.

044CAIN • 2 STAINLESS STEEL IMPRESSION COPINGS For OT CAP Normal and OT EQUATOR



Ref

044CAI22 • 2 STAINLESS STEEL IMPRESSION COPINGS Ø 2.25mm Spheres with interchangeable cap



Ref

• 2 STAINLESS STEEL IMPRESSION COPINGS
For OT CAP Micro

#### **TOOLS**



772CSF

HEX DRIVER - 0.9 mm For Threaded Micro Sphere

#### SINGLE THREADED SPHERES

#### WITH THREADED BONDING SLEEVE



Ref. 139KSFN Titanium + TiN Threaded Sphere With Sleeve For Bonding Kit - NORMAL SIZE

- 2 TITANIUM SINGLE THREADED SPHERES
- 1.3 mm Hex, 1.6 mm Thread 2 TITANIUM THREADED SLEEVES
- For Bonding 2 WAXING SPACERS

For Threaded Sphere - Normal Size

# Ref.: 139KSFM

Titanium + TiN Threaded Sphere With Sleeve For Bonding Kit - MICRO SIZE

#### Kit contains

- 2 TITANIUM SINGLE THREADED SPHERES
- 0.9 mm Hex, 1.6 mm Thread 2 TITANIUM THREADED SLEEVES
- For Bonding 2 WAXING SPACERS
- For Threaded Sphere Micro Size

#### SINGLE THREADED SPHERES

**NORMAL - MICRO** 

Ref.

039SFN2 •1 TITANIUM + TIN THREADED SPHERE NORMAL 1.3 mm Hex. 2.0 mm Thread



Ref :

Ref.:

Ref.:

039SFM2 •1 TITANIUM + TIN THREADED SPHERE MICRO 0.9 mm Hex, 2.0 mm Thread

#### **OT LOCK**



OT LOCK KIT

Kit contains: 880CLT 1 COMPLETE OT LOCK

- 1 BRASS POSITIONER
- 1 CERAMIC PIN



ADJUSTABLE OT LOCK KIT

#### Kit contains:

- 1 COMPLETE ADJUSTABLE OT LOCK • 1 EXTENDED BRASS POSITIONER
- 1 CERAMIC PIN
- 9 CASTABLE SPACER RINGS

#### OT LOCK SPARE PARTS



Ref

 CONICAL GUIDE 882CG



Ref.

882CAS • UNLOCKING TOOL

### INCLUDES OT CAP & OT BOX - OT STRATEGY - OT BAR - OT VERTICAL - OT UNILATERAL - OT EQUATOR

"BASIC" PROMOTIONAL KIT FOR LABORATORY

# Ref 005SKLBUS

#### TOOLS:

- 1 TWEEZER
- 1 PARALLELOMETER MANDREL OT CAP NORMO
- 1 PARALLELOMETER MANDREL OT CAP MICRO
- PARALLELOMETER MANDREL OT STRATEGY
- 1 PARALLELOMETER MANDRELOT BAR MULTIUSE
- 1 BLUE PLASTIC UNIVERSAL INSERTION HANDLE
- 1 INSERTION TOOL OT CAP NORMAL/MICRO
- **INSERTION TOOL OT STRATEGY**
- INSERTION TOOL OT BAR MULTIUSE
- INSERTION TOOL OT VERTICAL

#### Kit contains:

#### OT CAP - OT BOX:

- ASSORTED CASTABLE PIVOTS NORMAL / MICRO 16
- CASTABLE SPHERES NORMAL / MICRO • 4
- 2 CASTABLE OT CAP BARS NORMAL / MICRO
- CASTABLE OT BOX BARS CLASSIC (top + bottom) NORMAL / MICRO
- CASTABLE OT BOX SPECIAL BARS NORMAL / MICRO
- CASTABLE OT BOX CONNECTORS • 6
- CASTABLE OT BOX MONO HOUSING NORMAL / MICRO
- 8 POSITIONER RINGS NORMAL / MICRO
- 28 OT CAP RETENTIVE CAPS NORMAL / MICRO YELLOW, PINK, CLEAR, GREEN
- BLACK CAPS FOR PROCESSING NORMAL / MICRO q
- STAINLESS STEEL HOUSINGS NORMAL / MICRO FOR RESIN

#### OT STRATEGY:

- OT STRATEGY MALES 2 STANDARD BASE 2 LONG BASE
- CASTABLE STEADY

OT STRATEGY CAPS FOR STAINLESS STEEL **HOUSING** 

YELLOW, PINK, CLEAR

#### OT STRATEGY CONTINUED:

- 2 OT STRATEGY STAINLESS STEEL HOUSINGS
- 2 OT STRATEGY PLASTIC POSITIONING RINGS
- 6 OT STRATEGY CAPS FOR DUPLICATION **TECHNIQUE**

YELLOW, PINK, CLEAR

**BLACK CAPS - PROCESSING** (for wax and for duplication technique)

#### OT FQUATOR:

- CASTABLE MALES
- 2
- STAINLESS STEEL HOUSINGS RETENTIVE CAPS 2 PINK, 2 CLEAR
- BLACK CAPS FOR LABORATORY USE

#### OT UNILATERAL

- CASTABLE ATTACHMENT WITH COMBINED SPHERES
- **CASTABLE UNI-BOX**
- 1 MICRO POSITIONING RING
- 2 OT CAP MICRO CAPS - 1 PINK, 1 BLACK

- 2 OT STRATEGY CAPS - 1 PINK, 1 **BLACK For Duplication Technique**
- 1 CONNECTOR

#### OT BAR MULTIUSE:

- 1 CASTABLE BARS
- BAR EXTENSION
- POSITIONING CLIPS (Type A Type B)
- **CASTABLE BOXES** . 2
- CLIPS 2 PINK, 2 YELLOW • 4

#### OT VERTICAL:

- CASTABLE MALES • 2
- CASTABLE STEADY • 2
- 4 CLIPS - 2 WHITE, 2 GREEN PARALLELOMETER KEYS + PIN • 2
- CERAMIC PINS • 2



### KITS AND CODES

#### **IMPLANTOLOGY**

#### SPHERO FLEX - BLOCK SYSTEM TITANIUM ATTACHMENTS FOR OVERDENTURES

109

Ref.:

002



#### SPHERO FLEX

- 1 Titanium Abutment with self-aligning 2.5mm sphere 2 Pink Caps Soft Retention 1 Stainless Steel Housing
- 1 Protecitve Disk
- 3 Directional Rings

#### SPHERO BLOCK NORMAL

- 1 Titanium Abutment with stationary 2.5mm sphere 2 Pink Caps Soft Retention 1 Stainless Steel Housing
- Protective Disk
- 3 Directional Rings



#### SPHERO BLOCK MICRO

- 1 Titanium Abutment with stationary 1.8mm sphere 2 Pink Caps Soft Retention 1 Stainless Steel Housing 1 Protective Disk
- 3 Directional Rings

#### **ANCILLARY ITEMS**

Ref.:

003



Ref : 00PB MINI PARALLELOMETER

WITH UNIVERSAL TILTING MODEL TABLE (FOR LABORATORY USE, COURSES, ETC.)



**OT CEM** is a self and photo curing cement. It is designed for permanent metal to metal bonding in the use of attachments in prosthetic implant solutions

#### SPHERO FLEX / SPHERO BLOCK TOOLS



771CEF

UNIVERSAL KEY FOR SPHERO FLEX AND SPHERO BLOCK - NORMAL / MICRO



760CBM

Hex 23 mm



HEX DRIVER

FOR CONTRA-ANGLE TORQUE CONTROLLER



Ref. 760CBR

SCREW DRIVER FOR OT REVERSE Hex 1.3 mm

#### SPECIALTY ITEMS FOR IMPLANTS



### **CUFF HEIGHT MEASURING TOOL** Kit contains:

- 1 CUFF HEIGHT SLIDER GAGUE
  1 CUFF HEIGHT FIXED ROD GAGUE
  1 SILICON RINGS DISPENSER
- 20 SILICON RINGS



Ref.: 680FL 1 FCLAW REAMER BUR

#### **BROKEN SCREW EXTRACTOR KIT**

For removing broken screws from implants

contains

- 1 MANUAL CENTERING DEVICE 1 POSITIONER 1 CLAW REAMER BUR

- 1 REVERSE CUTTING BUR

Ref.: 680FS 1 REVERSE CUTTING BUR

Kit contains:

1 CLAW REAMER BUR
1 REVERSE CUTTING BUR

#### **IMPLANTOLOGY**

#### **ACCESSORIES FOR IMPLANTS**

For information on abutments for other implant systems please contact Rhein83



Ref : 108CV

Screw Vent Castable Abutment Non-Rotating with titanium screw White - Precision Hex 3.5 mm diameter



Ref. **108AVB** 

Screw Vent Castable Abutment Non-Rotating with titanium screw Red - Conical Hex For Bar Connections 3.5 mm diameter



108BRK

Branemark Castable Abutment Rotating with titanium screw 3.75 mm - 4.0 mm diameter



108BRK-NR

**Branemark Castable Abutment** Non-Rotating with titanium screw 3.75 mm - 4.0 mm diameter



108PE

Pitt Easy Castable Abutment Non-Rotating with titanium screw 3.25 mm - 3.75 mm - 4.0 mm diameter



108BFT

Straumann ITI Castable Abutment -Rotating with titanium screw for bar connections



Ref. 113BFT

Steel Transfer Abutment For Straumann ITI Implant with titanium screw



FA004

Steel Analog For Straumann ITI Implant

Rhein83 manufacutres castable abutments and titanium screws for most implant systems. For implant systems that are not listed in this catalog, please contact Rhein83 for additional information.



### RHEIN83 WORLD WIDE



#### A GLOBAL VISION WITH A COMMON TARGET

Our mission is to offer to the professionals of the dental field, different quality solutions allowing to reach the patient's comfort and satisfaction regardless the different social and financial situations.

This is possible to the precious support of our partners worldwide!



Rhein83 USA branch is active in the area since the year 2000 by supporting the distribution in the entire country including different areas in Latin America. Rhein83 USA is located in New Rochelle (few minutes away from NYC), taking care of developing an intense program of formation with courses dedicated to dentists and dental technicians. Courses will allow the attendants to have CTD's credits with speakers members of the "National Board for Certification in Dental Laboratory Technology, Inc":





Via ZAGO, 10/ABC 40128 - BOLOGNA

Tel. (+39) 051 244510 - (+39) 051 244396 Fax (+39) 051 245238

http://www.rhein83.com

Distributor





